

# SyMRI®

SyntheticMR AB

2019-12-16

## User Manual

Version 11.2.2

UDI: (01)07340024700030(8012)11.2.2

SyMRI 11

## Applicable Regulations

This section contains important information about regulatory clearance and usage of SyMRI. Note that SyMRI is classified as a medical device (MD) and is therefore subject to regulation that applies for clinical and non-clinical use. For instance, in the European union the Medical Device Directive (MDD) applies.

*SyMRI may only be used if and when legally permitted, check with local authorities and SyntheticMR AB if uncertain if this version of SyMRI can be legally used by you in your domain.*

---

**WARNING: If the regulatory domain that applies for you is not listed in this section SyMRI may not be used for any other purpose than demonstration or investigation as a non-medical device and shall not be used for diagnosis or any kind of patient management purposes.**

---



The marking “Medical Regulations Apply (Page I)” on each page of this manual is a reference to this section for further details.

### Certified Use

SyMRI complies with the following regulatory requirements. Usage is subject to local law and the conditions of the regulations. Additional restrictions, conditions or regulations may apply.

- CE (European Union):  
SyMRI 11 is a CE-marked medical device that complies with the Council Directive 93/42/EEC (MDD).
- USA/HSS (FDA): SyMRI 11 is an FDA 510(k) cleared medical device, K191036. Rx only.
- Australia (Therapeutic Goods Administration, TGA):  
For the Australian market the regulatory sponsor is:  
**Emergo Australia**  
201 Sussex Street, Level 20, Tower II Darling Park  
Sydney, NSW 2000  
Australia



### Not For Clinical Use (NFCU)

This section lists domains in which SyMRI may *under certain restrictions* be used in a non-clinical setting (NFCU). For domains where NFCU usage may apply, pay attention to the warning text that applies and any additional restrictions, conditions or regulations that may apply.

- Canada and the countries where Health Canada regulations is applicable:



**Investigational Device - To Be Used by Qualified Investigators Only**  
**Instrument de recherche - Réservé uniquement à l'usage de chercheurs compétents**

## Incident Reporting

In case of a serious incident, please contact the applicable representative, sponsor, and/or regulatory body that is applicable to the regulations where you have purchased and use SyMRI.

A complete list and additional legal information can be found at the following site:

<https://www.syntheticmr.com/legal>

## Patents and Trademarks

SyntheticMR AB's unique technology supports a faster MR imaging workflow and allows users to follow disease progression and therapy efficacy with greater confidence. A single 6-minute scanning sequence measures the absolute properties of the brain and delivers up to 8 different contrast weighted images, automatic segmentation and volumetric measurements of the patient. SyntheticMR AB is leading the way to a faster and more precise diagnosis of patients worldwide which has resulted in a patent portfolio covering many important innovations within synthetic MRI.

Please refer to the following page for an up-to-date listing of the applicable patent and trademarks, and additional legal information:

<https://www.syntheticmr.com/legal>

# Contents

Applicable Regulations .....	I
Certified Use .....	I
Not For Clinical Use (NFCU) .....	I
Incident Reporting .....	II
Patents and Trademarks .....	II
<b>Contents .....</b>	<b>II</b>
<b>1 Introduction .....</b>	<b>1</b>
1.1 Device Description .....	1
1.2 Intended Use .....	2
1.2.1 Indications for use .....	2
1.3 Advantages of SyMRI .....	3
1.4 Abbreviations .....	4
1.5 Warnings, Cautions and Notes .....	6
1.5.1 Note .....	6
1.5.2 CAUTION .....	6
1.5.3 WARNING .....	6
<b>Part I - Work with SyMRI® .....</b>	<b>6</b>
<b>2 Deployments of SyMRI .....</b>	<b>8</b>
2.1 SyMRI for Sectra PACS IDS7 (Plug-In) .....	9
2.1.1 Starting and Exiting .....	9
2.2 SyMRI StandAlone for Microsoft Windows .....	10
2.2.1 Starting and Exiting .....	10
2.3 SyMRI StandAlone for macOS .....	11
2.3.1 Starting and Exiting .....	11
2.3.2 Keyboard Shortcuts .....	11
2.4 SyMRI for Network .....	12
2.4.1 Starting and Exiting .....	12
2.5 SyMRI for syngo.via OpenApps .....	13
2.5.1 Starting and Exiting .....	13
<b>3 Concepts and Functions .....</b>	<b>14</b>
3.1 Concepts .....	14
3.1.1 Package .....	14
3.2 Workspace .....	15
3.2.1 Viewport .....	15
3.3 Toolbar .....	17
3.4 Options and Preferences .....	18
3.4.1 Layout Options .....	18
3.4.2 View Options .....	18
3.4.3 Other .....	19
3.5 Viewing Images .....	20
3.5.1 Zoom .....	20
3.5.2 Navigate .....	20
3.5.3 Windowing .....	20
<b>4 Optional Tools and Features .....</b>	<b>21</b>

4.1	AutoROI Tool	21
4.1.1	Functionality Description	22
4.2	Distance Measurement Tool	23
4.3	Quantification Table (Q-Table)	24
4.3.1	Displaying the Quantification Table	24
4.3.2	Information in the Quantification Table	24
<b>Part II - SyMRI® IMAGE</b>		<b>24</b>
<b>5 Functions in SyMRI IMAGE</b>		<b>26</b>
5.1	Synthetic Contrast Weighted Images	27
5.1.1	Display of Synthetic Contrast Weighted Images	27
5.1.2	Adjusting the Image Contrast	28
5.1.3	Adjusting FLAIR Images	29
5.1.4	Adjusting Double IR Images	30
5.1.5	Navigation Window for Contrast Images	30
5.2	SyMaps	32
5.2.1	Display of SyMaps	32
5.3	Region of Interest (ROI)	33
5.3.1	Create and Adjust ROI	33
5.3.2	Information in ROI	33
5.3.3	ROI Preferences	34
5.3.4	Functions on a ROI	34
5.4	Plots	36
5.5	Save and Load	37
5.5.1	Save SyMaps	37
5.5.2	Save a Screenshot of the Workspace	38
5.5.3	Save Image Stacks	38
<b>Part III - SyMRI® NEURO</b>		<b>39</b>
<b>6 Functions in SyMRI NEURO</b>		<b>41</b>
6.1	The Intracranial Mask	42
6.1.1	Adjusting the Intracranial Mask	42
6.1.2	Additional Functions for the Intracranial Mask	43
6.2	Tissue Segmentation	44
6.2.1	Brain Parenchyma and CSF Segmentation	44
6.2.2	Myelin Correlated (MyC) Partial Volume Mapping	44
6.3	Segmentation Maps	47
6.3.1	Display Segmentation Maps	47
6.3.2	Background Image for Segmentation	48
6.4	User Mask for User Defined Segmentation	49
6.4.1	Display the User Mask	49
6.4.2	Measuring Volume Using the User Mask	49
6.4.3	Create and Adjust User Mask	49
6.4.4	Save and Load the User Mask	50
6.5	Segmentation Table	52
6.5.1	Displaying Volumes per Slice	52
6.5.2	User Mask in the Segmentation Table	53
6.5.3	Exporting the Segmentation Table	53
6.6	Reference curves	54

6.7 ROI in SyMRI <i>NEURO</i> .....	55
6.7.1 Measuring the Volume of Segmented Tissue .....	55
6.8 Plots in SyMRI <i>NEURO</i> .....	56
<b>Part IV - Important Information .....</b>	<b>56</b>
<b>7 Warnings .....</b>	<b>58</b>
7.1 Image Artifacts .....	58
7.1.1 Motion Artifacts .....	58
7.2 Deviations - Synthetic MRI .....	59
7.2.1 Black Blood .....	59
7.2.2 Flow Sensitivity .....	59
7.2.3 Small Features Might Disappear .....	59
7.3 Functions Specific for Synthetic MRI .....	61
7.3.1 Use ROI to Null Tissue .....	61
7.3.2 Quantitative Measurements .....	61
7.3.3 Segmentation .....	62
7.3.4 Factors Affecting the Volume Estimations .....	63
7.4 Additional Warnings .....	65
7.4.1 System Requirements .....	65
7.4.2 Follow-up Examinations Using Different Analysis Tools ..	66
7.4.3 Character Encoding .....	66
7.4.4 Anonymization Software .....	66
7.5 Scanner Specific Warnings .....	67
7.5.1 GE .....	67
7.5.2 Philips .....	67
7.5.3 Siemens .....	68
<b>Part V - Appendix .....</b>	<b>68</b>
<b>8 Keyboard Shortcuts .....</b>	<b>70</b>

# 1 Introduction

## 1.1 Device Description

SyMRI allows the user to generate multiple image contrasts from a single acquisition scan. This is accomplished by post-processing a multi-delay, multi-echo acquisition (MDME) into parametric maps. The parametric maps are R1, R2 relaxation rates, and proton density (PD). The inverse relaxation parameters, T1 relaxation time ( $1/R1$ ), and T2 relaxation time ( $1/R2$ ) are also provided as parametric maps. The parametric maps can be visualized as contrast weighted MR images, such as T1-, T2-, PD, and Inversion Recovery (IR) weighted images (including T1W FLAIR, T2W FLAIR, STIR, Double IR, and PSIR weighted images). SyMRI calculates the pixel signal intensity as a function of R1, R2, PD, and desired MR scanner settings, such as echo time (TE), repetition time (TR), and inversion delay time (TI). A number of default settings for TE, TR, and TI are provided, but the user has the ability to change the contrast of the images. SyMRI generates all the different image contrasts from the same parametric maps, derived from the same acquisition. This leads to enhanced image slice registration, owing to the absence of inter-acquisition patient movement. SyMRI provides the user the ability to change the contrast of the images after the acquisition. This is performed by adjusting the TE, TR, and/or TI parameters post-acquisition, to generate the specific contrast desired.

SyMRI also enables the users to obtain volumetric information in the head, including white matter (WM), grey matter (GM), cerebrospinal fluid (CSF), Myelin correlated (MyC) partial volume, brain parenchyma (BP) and intracranial cavity (IC). This is accomplished by using tissue definitions based on the parametric maps. The tissue definitions provide tissue partial volume, or tissue fraction, per voxel. SyMRI also provides image processing tools to extract the values of the parametric maps, and tissue partial volume, per individual pixel, per region of interest, or the entire imaging volume.

SyMRI is intended to be used on data produced by any of the following acquisition sequences:

- MDME sequence data from GE MAGiC
- MDME sequence data from Philips SyntAc
- MDME sequence data from Siemens 3 T TSE\_MDME
- MDME sequence data from Siemens 1.5 T TSE\_MDME<sup>1</sup>

---

<sup>1</sup> Not yet 510(k) cleared by the US FDA and thus currently not available for sale in the US unless limited to research or investigational use

## 1.2 Intended Use

### 1.2.1 Indications for use

SyMRI is a post-processing software medical device intended for use in visualization of the brain. SyMRI analyzes input data from MR imaging systems. SyMRI utilizes data from a multi-delay, multi-echo acquisition (MDME) to generate parametric maps of R1, R2 relaxation rates, and proton density (PD). SyMRI can generate multiple image contrasts from the parametric maps. SyMRI enables post-acquisition image contrast adjustment. SyMRI is indicated for head imaging.

SyMRI is also intended for automatic labeling, visualization and volumetric quantification of segmentable brain tissues from a set of MR images. Brain tissue volumes are determined based on modeling of parametric maps from MDME images.

When interpreted by a trained physician, SyMRI images can provide information useful in determining diagnosis. SyMRI should always be used in combination with at least one other, conventional MR acquisition (e.g. T2W FLAIR).



## 1.3 Advantages of SyMRI

SyMRI only uses a single acquisition to obtain R1, R2 and PD values. Synthetic T1W, T2W and FLAIR images with any choice of echo time TE, repetition time TR and inversion delay time TI can be shown without the need to rescan the patient. Therefore, the method can save valuable scan time for the MR examination.

**CAUTION:** The synthetic T2W FLAIR shows an edge-enhancement effect at the interface of CSF and brain tissue. It is the result of partial volume effects of the acquisition voxels at the interface.

Use of the synthetic T2W FLAIR contrast depends on the clinical question. The synthetic T2W FLAIR has an excellent sensitivity for changes in brain tissue.

**CAUTION:** We recommend that any changes visible in T2W FLAIR are confirmed using the T2W, PDW, PSIR or Double IR of the same synthetic dataset.

Since all synthetic contrast weighted images are generated from the same parametric maps, from the same acquisition, the synthetic images will be registered. This means that a finding in T2W FLAIR will appear in the exact same pixel in the other synthetic contrast weightings (T2W, PDW, PSIR or Double IR) for confirmation or rejection as a false positive.

**CAUTION:** For diagnosis of illnesses resulting in edge-enhancements in the T2W FLAIR, such as inflammatory reactions or hemorrhages, the synthetic sequence alone may not be sufficient for diagnosis. In case of doubt, the addition of for example a conventional 3D-FLAIR is recommended.

---

**WARNING:** Patient management decisions should not be made based solely on synthetic contrast weighted images from SyMRI and should always be used in combination with conventional MRI or other diagnostic tools applicable.

---



The tissue segmentation features provide quantitative measures. This can support objective assessment of the patient. Multi-slice imaging and enhanced processing techniques using the SyMRI software generates a three-dimensional volumetric quantification of brain, CSF, WM, GM, MyC and total intracranial volumes.

This user manual provides an overview of available functionality and describes all concepts used in the software.

## 1.4 Abbreviations

The following abbreviations are used in the manual.

BPV	The Brain Parenchymal Volume; the volume of all tissues in the ICV excluding the CSF volume
BPF	The Brain Parenchymal Fraction, the BPV divided by the ICV
CSF	SyMRI's definition of Cerebrospinal Fluid (see Section 6.2)
DIR	Double IR, an inversion recovery image with two inversion times
FLAIR	FLuid Attenuated Inversion Recovery, an MR image where the signal from CSF is suppressed
GM	SyMRI's definition of gray matter (see Section 6.2)
HUD	Heads-Up Display; the information presented alongside an image in SyMRI, like patient- and image information
ICV	Intracranial Volume; the total volume of the Intracranial mask
IR	Inversion Recovery; imaging where an inversion prepulse is applied.
MRI	Magnetic Resonance Imaging
MyC	SyMRI's tissue definition that correlates with Myelin
Non-WM/GM/CSF	Tissue with characteristics which deviate from WM, GM and CSF (see Section 6.2)
PD	The proton density
PDW	An image with PD-weighted contrast, typically short TE and long TR
PSIR	Phase Sensitive Inversion Recovery, an image with inversion recovery where the signal sign is kept
R1	The longitudinal R1 relaxation rate (1/T1)
R2	The transverse R2 relaxation rate (1/T2)
ROI	Region of Interest; a user-defined region in one slice in the data set, drawn with a dedicated ROI-tool
STIR	Short T1 IR image, a fat suppression technique
SyMaps	A map that shows one of the quantitative values T1, T2, R1, R2 or PD
T1	The longitudinal T1 relaxation time (1/R1)
T1W	An image with T1-weighted contrast, typically short TE and short TR
T2	The transverse T2 relaxation time (1/R2)
T2W	An image with T2-weighted contrast, typically long TE and long TR
TE	The echo time
TI	The inversion delay time after an inversion pulse

TR	The repetition time
UM	User Mask, see Section 6.4
WM	SyMRI's definition of (myelinated) white matter (see Section 6.2)

## 1.5 Warnings, Cautions and Notes

### 1.5.1 Note

Notes contain important information about features in SyMRI that needs to be highlighted. Information presented as Notes includes helpful tips, specific limitations to the availability of certain functions, as well as highlighting behavior in SyMRI that may not be intuitive to the user.

Note: Example of a note

### 1.5.2 CAUTION

Cautions indicate potentially hazardous situations for the patient or user that could result in lost time, reduced image quality and/or the patient having to be re-examined. Cautions also indicate how such situations are avoided.

CAUTION: Example of a caution

### 1.5.3 WARNING

Warnings indicate potentially hazardous situations that could result in injury for the patient, mainly in the form of misinterpretation and misdiagnosis.

Warnings also indicate how such situations are avoided.

The ISO 15223-1:2016 warning symbol is used with these warnings to highlight the importance of the warning.

---

WARNING: Example of a warning

---



# Part I

Work with SyMRI®

## 2 Deployments of SyMRI

SyMRI is available in several deployments, depending on your system and licensing one or more deployment modes are available. Please consult your system administrator or manager for queries on which deployments are available for your use.

Most functionality, behavior and user interface interactions are identical across all deployments whereas some aspects differ slightly. The following chapters of this part of the user manual describes the differences and some particulars of the deployments.

In some places in the other parts of the user manual there are sections that only apply to one of the deployments. These sections are marked with a margin note and a vertical bar in the following manner:

**Plug-In**

This paragraph only applies to SyMRI for Sectra PACS IDS7 (Plug-In).

**StandAlone**

This paragraph only applies to SyMRI StandAlone.

**Windows**

This paragraph only applies to SyMRI for Microsoft Windows.

**macOS**

This paragraph only applies to SyMRI for macOS.

**Network**

This paragraph only applies to SyMRI for Network.

## 2.1 SyMRI for Sectra PACS IDS7 (Plug-In)

When SyMRI is used as a plugin in Sectra PACS IDS7 the application toolbar and menu are placed in a floating toolbar window instead of attached to the main window. The floating toolbar can be freely moved at your discretion.

The workspace with the floating toolbar can be seen in Figure 2.1.

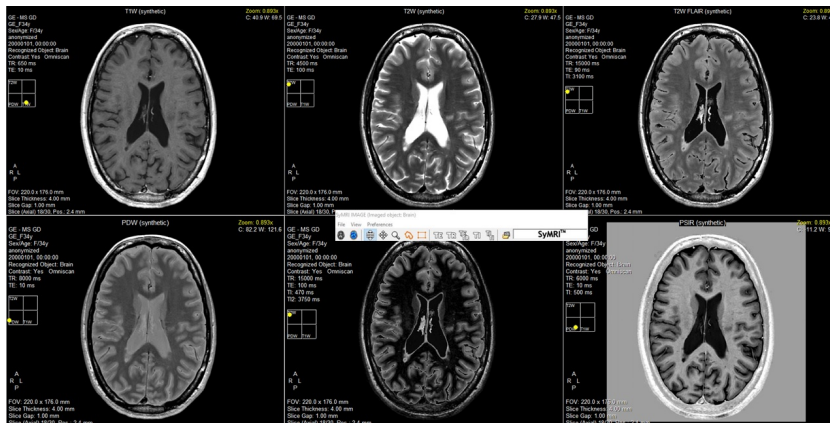


Figure 2.1: The workspace in SyMRI with a floating toolbar.

### 2.1.1 Starting and Exiting

This chapter describes how to start SyMRI on a Sectra IDS7 workstation.

#### 2.1.1.1 Start SyMRI

- Start and log on to Sectra IDS7.
- Open desired image stack in a Sectra IDS7 image window.
- Select “SyMRI” in the clinical application menu.
- SyMRI starts and can be used.

In the event of non-recoverable errors during loading, SyMRI will terminate and return to Sectra PACS IDS7 after displaying a message box with information about the error.

#### 2.1.1.2 Close SyMRI

SyMRI is closed by using the menu item File>Exit SyMRI in the floating toolbar (see Section 3.3).

Do not close the Sectra image window before exiting SyMRI since it can cause incorrect termination of SyMRI.

## 2.2 SyMRI StandAlone for Microsoft Windows

### 2.2.1 Starting and Exiting

This chapter describes how to start the StandAlone deployment of SyMRI.

#### 2.2.1.1 Start SyMRI

- Start the SyMRI stand-alone program.
- Click on the main window background and select the desired dataset.
- SyMRI starts and can be used.

In the event of a dataset related error during loading, SyMRI will return to the main window where you can click to start on a different dataset after a message box with information about the error has been displayed. If the error is not recoverable, SyMRI may exit after displaying the error message.

#### 2.2.1.2 Close SyMRI

SyMRI can be closed using the menu item File>Exit SyMRI in the toolbar (see Section 3.3) or by pressing the main window's close button.



## 2.3 SyMRI StandAlone for macOS

### 2.3.1 Starting and Exiting

This chapter describes how to start the StandAlone deployment of SyMRI.

#### 2.3.1.1 Start SyMRI

- Start the SyMRI stand-alone program.
- Click on the main window background and select the desired dataset.
- SyMRI starts and can be used.

In the event of a dataset related error during loading, SyMRI will return to the main window where you can click to start on a different dataset after a message box with information about the error has been displayed. If the error is not recoverable, SyMRI may exit after displaying the error message.

#### 2.3.1.2 Close SyMRI

SyMRI can be closed using the menu item SyMRI>Exit SyMRI in the menubar (see Section 3.3) or by pressing the main window's close button.

### 2.3.2 Keyboard Shortcuts

On macOS the following shortcut keys are different from Windows:

- ⌘ instead of Ctrl.
- F13 instead of Insert.

## 2.4 SyMRI for Network

### 2.4.1 Starting and Exiting

This chapter describes how SyMRI is used in networking mode.

#### 2.4.1.1 Start SyMRI

- Start SyMRI from PACS.
- Depending on if a single series or a study has been selected, a list of series which can be loaded is presented.
- If there are saved user masks or intracranial masks, one of these may be selected after the series has been loaded from the network.

In the event of non-recoverable errors during loading, SyMRI will terminate after displaying a message box with information about the error.

#### 2.4.1.2 Close SyMRI

SyMRI can be closed using the menu item File>Exit SyMRI in the toolbar (see Section 3.3) or by pressing the main window's close button.

## 2.5 SyMRI for syngo.via OpenApps

### 2.5.1 Starting and Exiting

This chapter describes how SyMRI is used in syngo.via OpenApps.

#### 2.5.1.1 Start SyMRI

- Start SyMRI from syngo.via OpenApps. See the User Manual for syngo.via OpenApps for more information.
- If a series has previously been opened with SyMRI, double click on the series for syngo.via OpenApps to automatically start SyMRI with the series.

In the event of non-recoverable errors during loading, SyMRI will terminate after displaying a message box with information about the error.

#### 2.5.1.2 Close SyMRI

SyMRI can be closed by pressing the tab's close button. See the User Manual for syngo.via OpenApps for more information.

## 3 Concepts and Functions

This chapter describes general concepts, functions, user commands and functions accessible in all packages within SyMRI. Specific functionality for each separate package is covered in the following chapters.

### 3.1 Concepts

#### 3.1.1 Package

There are two packages available in SyMRI, SyMRI *IMAGE* (Chapter 5) and SyMRI *NEURO* (Chapter 6). Depending on your licensing one or more of them are available for use.

When a dataset has been loaded, the active package can (if your licensing permits) be changed via the menu item Preferences > Switch Package.

##### 3.1.1.1 SyMRI *IMAGE*

The SyMRI *IMAGE* package is available with an *IMAGE* or *NEURO* license. See Chapter 5 for a description of the SyMRI *IMAGE* package and its functionality.

##### 3.1.1.2 SyMRI *NEURO*

The SyMRI *NEURO* is only available with a *NEURO* license and will automatically be selected when a brain dataset is loaded and classified as a brain (the text: "Recognized Object: Brain" appears on the left side of the image).

If the dataset is classified as another anatomy/object, the SyMRI *IMAGE* package will be selected instead. You can change to SyMRI *NEURO* using the menu item Preferences > Switch Package.

SyMRI *NEURO* contains all features of SyMRI *IMAGE*, but also provides automatic tissue segmentation, volume calculations and BPF-measurement. See Chapter 6 for a description of the SyMRI *NEURO* package and its functionality.

## 3.2 Workspace

The complete window covered by SyMRI is called the workspace. The workspace (Figure 3.1) includes one or several viewports with images or tables.

The workspace is controlled by using the SyMRI menu and toolbar (section 3.3). The keyboard and mouse can also be used.

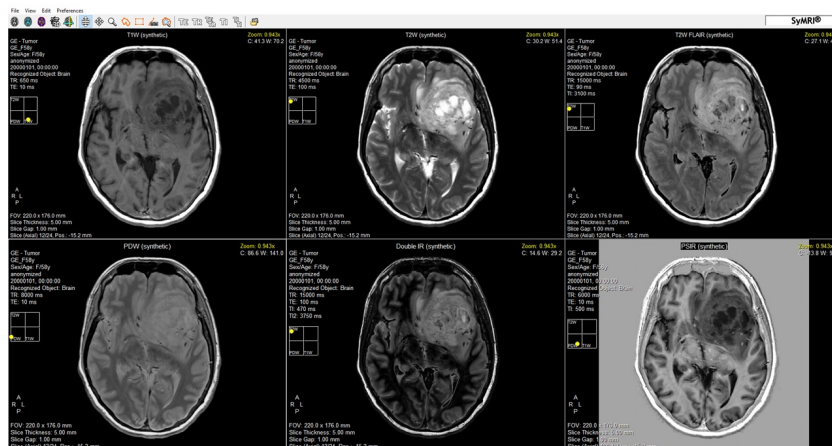


Figure 3.1: The workspace in SyMRI configured with six viewports. The contents in the toolbar and menus depend on the current package and your licensing.

### 3.2.1 Viewport

Each separate image or table within the workspace of SyMRI is called a viewport. Each viewport can display one image with image information (Figure 3.2), or (if your licensing permits) a table of tissue volumes. The viewport where the mouse pointer is placed is always active. The viewport content is controlled by the right-click menu.

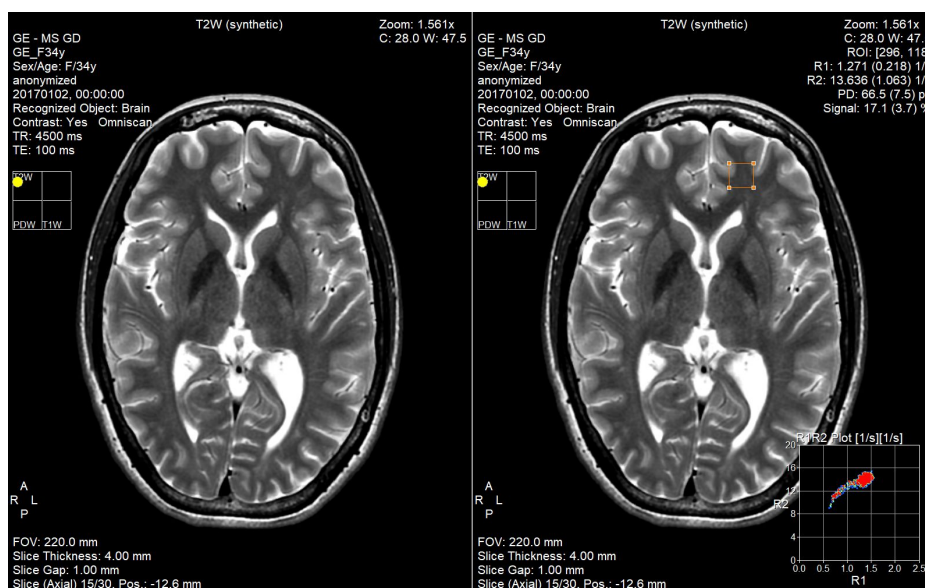


Figure 3.2: Two viewports in the SyMRI workspace. When a ROI is active, additional information is displayed on the right side of the viewport.

### 3.2.1.1 Viewport Image Information

The viewport image information can also be referred to as the “HUD” (Heads-Up Display). The information displayed in the HUD depends on what type of image is shown in the viewport.

**Patient information:** The patient information is displayed in the upper left corner and includes patient name, ID, age, and gender.

**Examination information:** The examination information is displayed in the upper left corner, below the *patient information*. Protocol name, acquisition date and time is displayed, and if applicable, information about administered contrast agent.

Note: If any of the displayed information is missing from the data set, e.g. due to anonymization, the line usually displaying this data will be completely removed from the HUD and the subsequent lines will be moved up so that no empty space is visible.

---

**WARNING:** Information about the contrast agent/bolus is shown *only if* the dataset contains the DICOM tag that indicates that a contrast agent or bolus was used in the scanning. When this tag is present, a line starting with “Contrast:” is shown. If the tag is present but has no value, the line will be “Contrast: Yes”. If the tag has a custom value, the custom value will be displayed instead of Yes. When this DICOM tag is not present, meaning that no contrast agent/bolus was used, the contrast/bolus line is not shown.

---



**Acquisition geometry:** The acquisition geometry is displayed in the bottom left corner of the viewport. This includes FOV, slice thickness, slice gap, number of slices and the position of the current slice.

**Zoom and Windowing:** The current zoom, as well as the center and width of the windowing is displayed in the upper right corner.

Information on contrast weighting is only visible in viewports where a contrast weighted image is displayed.

**Contrast weighting information:** The contrast weighting information is displayed in the upper left corner, below the *Examination information*, and at the top of the viewport. The information to the left contains the scanner parameters, like TE, TR, etc. (Section 5.1). The current contrast weighting is written at the top of the viewport.

**Navigation window:** The navigation window appears below the text in the upper left corner. (Section 5.1.5)

There is also some information that is only displayed when a ROI is active in a viewport (Section 5.3).

**ROI info:** Statistical information about the ROI appear in the upper right corner, below the *Zoom and Windowing* information. It contains information about the tissue properties within the ROI, as well as the signal intensity (Section 5.3.2).

**ROI plot:** A scatter plot also appears in the bottom right corner of the viewport (Section 5.4).

Some additional information may be displayed when using a SyMRI *NEURO* license (Section 6).

### 3.3 Toolbar

The commands available in the SyMRI toolbar target the entire workspace rather than one specific viewport. The contents in the toolbar and menus depend on the current package and your licensing. For example, language and workspace settings are found in the toolbar, along with functions for exporting entire datasets and screenshots.

The availability of some features depends on your licensing.



Figure 3.3: The toolbar in SyMRI. The upper is the toolbar for SyMRI *NEURO* with the AutoROI Tool. The lower is the toolbar in SyMRI *IMAGE*.

The specific commands in the toolbar and on the right-click menu are described in the following sections.

#### Plug-In

Note: Even though SyMRI is integrated in Sectra IDS7, it is a separate program. The tools available in the Sectra IDS7 application *can NOT* be used in SyMRI and vice versa. Only the tools on the SyMRI floating toolbar are applicable to SyMRI.

Pressing the close button or closing the Sectra IDS7 main window without exiting SyMRI may cause SyMRI to terminate incorrectly or leave the floating toolbar window accessible until a new patient or examination is selected.

#### macOS

In macOS the menubar is found at the top of the screen when SyMRI is in focus.

## 3.4 Options and Preferences

### 3.4.1 Layout Options

The number of viewports can be changed.

**View>1 Viewport** «Ctrl+1» Show one viewport

**View>2 Viewports** «Ctrl+2» Show two viewports

**View>3 Viewports** «Ctrl+3» Show three viewports

**View>4 Viewports** «Ctrl+4» Show four viewports

**View>5 Viewports** «Ctrl+5» Show five viewports

**View>6 Viewports** «Ctrl+6» Show six viewports

Viewport layouts, including the settings of each viewport, can be saved and loaded.

**Preferences>Save Current Layout** Save current layout with number of viewports, time parameters, overlay type, type of SyMaps and image settings. When saving the layout, a dialog box asking for a name of the layout is shown.

If there are any saved layouts, SyMRI will select one as start-up layout whenever SyMRI is started. Which layout is selected depends on three criteria:

1. Acquisition orientation
2. Recognized object
3. Field strength

If no exact match can be found, the closest match will be selected by matching the criterias in the order specified above.

**Preferences>Load Layout** Load a saved layout. Each saved layout is an item in a submenu with the name given when saving the layout. Selecting one of the items loads that layout. This layout will be used as default layout when starting SyMRI.

**Preferences>Delete Layout** Delete a saved layout. Each saved layout is an item in a submenu with the name given when saving the layout. Selecting one of the items deletes that layout.

### 3.4.2 View Options

**View>Link Zoom and Pan** All viewports in the workspace are linked. The pan and zoom will be adjusted simultaneously in all viewports. This will make the zoom and pan identical for all images.

**View>Bicubic interpolation** «Ctrl+P» Interpolation gives a smoother visualization of contrast weighted images and SyMaps. Interpolation of images is active by default and is always linked between all viewports. This setting is not affected by, and does not effect, *View>Link Zoom and Pan*.

Note: Zooming out far while interpolation is turned off may cause aliasing artifacts in the image.

**Right-Click>Image Settings>Color** «C» When the color option is activated, the image stack is displayed in color instead of grayscale. Synthetic contrast weighted images cannot be visualized in color. (See Section 6.3.2.)



### 3.4.3 Other

**Preferences > Language** Choose desired language from the list of available languages. When any language other than English is used the menu item always contains the text "(Language)" (in English) after the word for language in the currently chosen language. The sub-menu contents are the names of each available language, written in that language.

#### Windows

**File > About SyMRI** Information about the currently running version of SyMRI will be displayed.

#### macOS

**SyMRI > About SyMRI** Information about the currently running version of SyMRI will be displayed.

## 3.5 Viewing Images

### 3.5.1 Zoom

**Zoom Tool** «Z» Moving the mouse while keeping the left mouse button down will zoom the image stack in the viewport. Movement forward zooms in and movement backward zooms out. Tables and plots cannot be zoomed.



**Right-Click>Zoom>Zoom in** «Ctrl+Z» The image is zoomed in one step.

**Right-Click>Zoom>Zoom out** «Ctrl+Shift+Z» The image is zoomed out one step.

The zoom steps are predefined.

**Right-Click>Zoom>Zoom to Fit Window** The size of the image is adjusted to the size of the viewport.

**Right-Click>Zoom>Pixel to Pixel** One voxel in the original DICOM-dataset is displayed using one pixel on the computer screen. For non-square voxels the shortest side of the voxel in the image plane is mapped to one pixel on the computer screen.

### 3.5.2 Navigate

**Scroll through image stack:** Scroll through the image stack by scrolling the mouse wheel or using the «Page Up» and «Page Down» keys. The first slice of an image stack is reached through «Home». The last slice is reached through «End».

The following tools can be found in the SyMRI toolbar.

**Paging Tool** «F4» Moving the mouse while keeping the left mouse button down will scroll through the image stack.



**Pan Tool** «F5» «P» Moving the mouse while keeping the left mouse button down will pan the image stack in the viewport. Tables and plots cannot be panned.



### 3.5.3 Windowing

**Windowing:** Keep the middle mouse button or mouse wheel down and move the mouse to adjust the contrast and intensity of the image. *Horizontal movement* will change the window centre and *vertical movement* changes the window width. The centre and width will determine the gray- or color-scale of the image, which in turn affects brightness and contrast.

**Right-Click>Image Settings>Auto Scale** «A» The image brightness and contrast is set to suit the image in the viewport.

**Right-Click>Image Settings>Invert Grayscale** «Ctrl+Shift+I» When the inverted grayscale option is activated, the image stack is displayed with the intensity of each pixel inverted i.e. black pixels are displayed as white and vice versa.

## 4 Optional Tools and Features

The availability of these optional tools may depend on your choice of licensing, please contact your sales representative for inquiries about purchasing optional features.

### 4.1 AutoROI Tool

A tool that combines and enhances functions of the ROI and user mask for quicker manual segmentation.

The ROI and user mask created by AutoROI Tool has the same properties and basic functionality as those created with the basic ROI and user mask functionality. See sections 5.3 and 6.7 for more details on basic ROI operations and features, and section 6.4 for the User Mask.

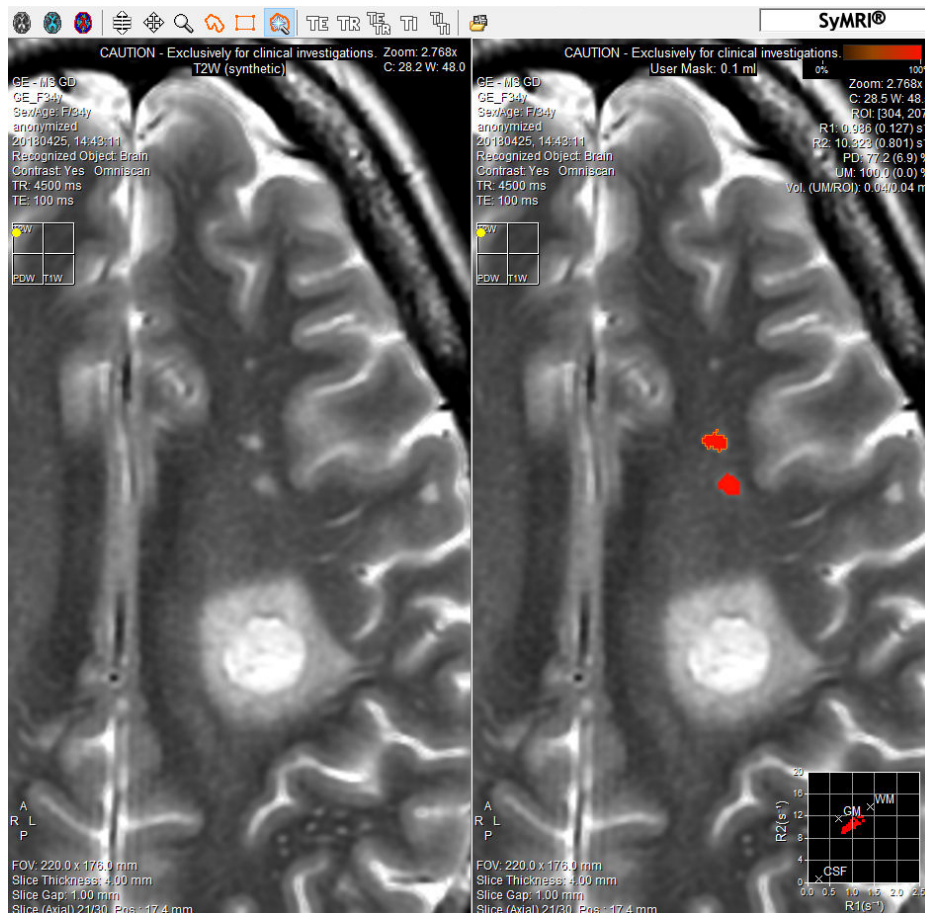


Figure 4.1: A screenshot of lesions segmented using the AutoROI Tool.

Note: The AutoROI Tool does not make any pathological classification of the selected features. Clinical competency is required for correctly assessing the characteristics of the ROI and adjusting the generated user mask if needed.

### 4.1.1 Functionality Description

The AutoROI Tool is intended to facilitate and make manual segmentation of features and structures in the brain easier and less time consuming.



**AutoROI Tool** «X» Click or draw to add or remove regions to or from the user mask.

**Toggle add/remove mode:** «X» When the AutoROI Tool is active, toggle if the AutoROI Tool should add to or remove from the user mask.

To activate the tool, click the AutoROI Tool icon in the toolbar or press «X» once. Pressing «X» after the tool has been activated will toggle between adding to the user mask and subtracting from the user mask. The tool is initially in the adding mode when enabled, when in the subtractive mode the mouse cursor changes with an minus sign indicating the subtractive mode (also see table 4.1).



Table 4.1: AutoROI Tool cursor for the two modes of operation.

Cursor Icon	AutoROI Mode
	Adding mode
	Subtracting mode

#### 4.1.1.1 Draw Functionality

A ROI of free shape can be drawn in the same way as a freehand ROI (section 5.3.1). When the ROI is drawn on a contrast weighted image, the contained pixels are automatically added to the user mask.

When the ROI is drawn on a segmentation overlay, the ROI will contain the selection within the overlay map, thus selecting only within the values segmented as the selected tissue/fluid.

In the addition mode the region is added to the user mask, and conversely in the subtractive mode, the region is subtracted from the user mask.

#### 4.1.1.2 Click Functionality

When the AutoROI Tool is used to click on features on contrast weighted images or segmentation overlay, the feature is automatically segmented and added to the user mask. A ROI is created around the segmented region. If the segmented region overlaps with existing user mask the ROI is created around the whole region.

In the addition mode the region is added to the user mask, and conversely in the subtractive mode, the region is subtracted from the user mask.

Note that when used on contrast weighed images, the functionality is primarily designed for easier selecting findings such as MS-lesions<sup>1</sup> and other small morphological changes. For larger morphological changes additional clicks may be required to sufficiently cover the entire morphology of interest.

<sup>1</sup> Additional manual adjustments of the selection might be required.

When used on a segmentation overlay, the function is primarily designed and optimized for ventricle segmentation. The tool can be used on any segmentation overlay available but might have less optimal behavior.

## 4.2 Distance Measurement Tool

The Distance Measurement Tool is a tool to measure the distance between the center of two voxels within a slice. It is activated by selecting the icon in the toolbar.



A distance is measured by placing the mouse pointer on the start point in the viewport, holding down the left mouse button and moving the pointer to the endpoint, and finally releasing the left mouse button. The distance is shown as a line with two endpoints and the distance in mm printed. The distance can be moved afterwards by holding down the left mouse button with the mouse pointer on the line and moving it. The distance can also be changed by moving only one of its endpoints by holding down the left button with the mouse pointer over one of the endpoints and moving it.

Right clicking on a distance activates a context menu:

**Right-Click>Disable Distance** Remove the distance clicked on.

**Right-Click>Disable All Distances** Remove all distances in the current viewport.

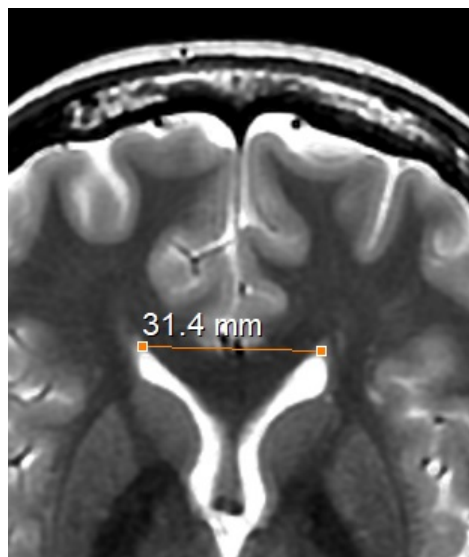


Figure 4.2: An example of using Distance Measurement Tool.

## 4.3 Quantification Table (Q-Table)

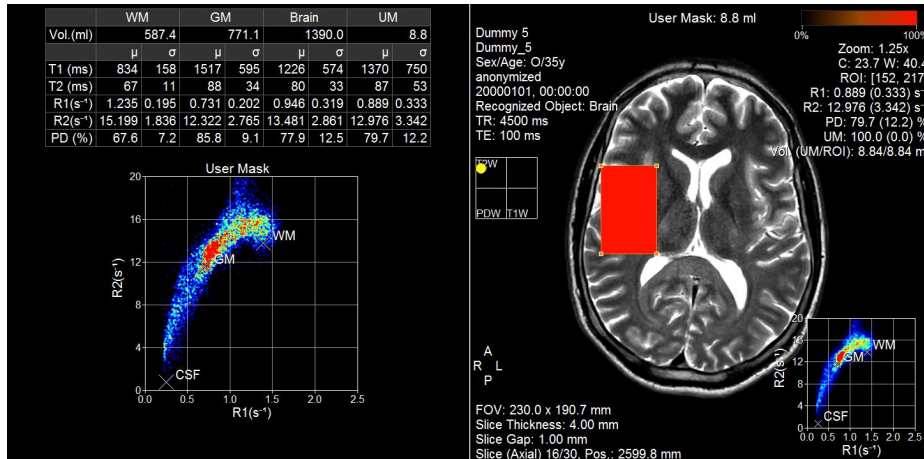


Figure 4.3: An example of the Quantification Table on a dataset with a user mask.

### 4.3.1 Displaying the Quantification Table

To display the Quantification Table, right-click and select “Quantification Table”.

### 4.3.2 Information in the Quantification Table

The Quantification Table is intended to give additional information about the quantitative values in WM, GM, Brain and User Mask. “Brain” is the brain parenchyma, the combination of WM, GM, and NON. Values displayed are mean and standard deviation of T1, T2, R1, R2, and PD within the segmented tissues. The total volumes are also displayed. Unlike the Segmentation Table, the Quantification Table doesn’t display individual slices. The user mask can be used to display the tissue values in chosen regions of one or several slices. See section 6.4 (User Mask for User Defined Segmentation) in the user manual for information about how to use the user mask. This section also describes how to copy a segmented tissue, such as CSF or myelin, to the user mask and thereby adding an additional tissue to the Quantification Table.

A larger version of the Q-plot is also displayed. The axis of the Q-plots can be changed to R1-PD and R2-PD in the same way as the regular Q-plot.

# Part II

SyMRI<sup>®</sup> *IMAGE*

## 5 Functions in SyMRI *IMAGE*

In SyMRI *IMAGE* the functionality to perform Synthetic Magnetic Resonance Imaging is available. The approach uses a single MR quantification scan to measure the tissue properties R1 and R2 relaxation rates and proton density PD. Based on these three properties contrast weighted images such as T1W, T2W or FLAIR can be synthesized with a free choice of echo time TE, repetition time TR, in case of an inversion pulse, inversion delay time TI, and in case of a double inversion pulse, a second inversion delay time TI2. Default settings for common contrast weighted images are provided but the user may alter TE, TR, TI and TI2 to any possible value.

**CAUTION:** The synthetic images are acquired and calculated differently than actual conventional images and therefore the image contrast is very similar, but not necessarily identical, even at the same settings for TE, TR and TI. Take your time to find a setting of TE, TR and TI such that either the appearance of the images is similar to those you are accustomed to, or optimized according to your needs. Interpretation of synthetic images should be done with care and an educational period to learn the nuance differences between synthetic and conventional images is recommended.

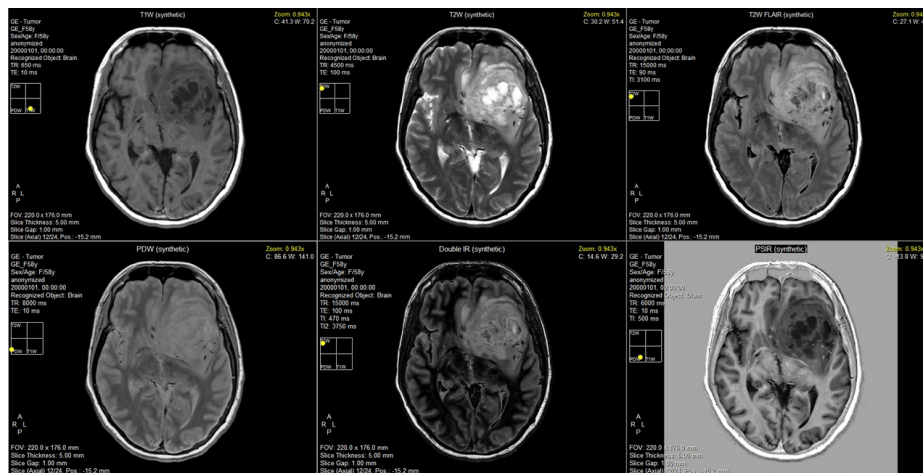


Figure 5.1: When SyMRI *IMAGE* is started, six viewports are shown with synthetic contrast weighted images: T1W, T2W, T2W FLAIR, PDW, Double IR (WM supp), and PSIR



## 5.1 Synthetic Contrast Weighted Images

The synthetic contrast weighted images are calculated using the measured R1, R2 and PD values in combination with the selected TE, TR, T1 and T12 settings.

### 5.1.1 Display of Synthetic Contrast Weighted Images

A number of preset contrast weightings can be selected using the right-click menu or keyboard commands. The available settings may depend on your licensing.

**Contrast Layout** «Ctrl+R» Six viewports with contrast images are displayed: a T1W, a T2W a T2W FLAIR, a PDW, a Double IR (WM supp) and a PSIR. Default values for TE, TR and TI are listed in table 5.1. Table 5.2 shows the defaults for DIR when a DIR enabled license is used.

**Right-Click>T1W** «1» A T1-weighted image is displayed in the active viewport. In T1-weighted images, areas with low proton density are weighted to enhance the image contrast between gray and white matter.

**Right-Click>T2W** «2» A T2-weighted image is displayed in the active viewport.

**Right-Click>PDW** A PD-weighted image is displayed in the active viewport.

**Right-Click>T1W FLAIR** A T1-weighted Fluid Attenuated IR image is displayed in the active viewport.

**Right-Click>T2W FLAIR** «3» A T2-weighted Fluid Attenuated IR image is displayed in the active viewport.

**Right-Click>STIR** «4» A T2-weighted Short Tau IR image is displayed in the active viewport.

**Right-Click>PSIR** «5» A Phase Sensitive IR image is displayed in the active viewport.

**Right-Click>PSIR (vessel)** A Phase Sensitive IR image to highlight blood vessels is displayed in the active viewport.

**Right-Click>Double IR (WM supp)** A Double IR image is displayed in the active viewport, which suppresses both WM and CSF.

**Right-Click>Double IR (GM supp)** A Double IR image is displayed in the active viewport, which suppresses both GM and CSF.

Note: The DIR menu items are only available if you have a license that enables the DIR feature.



Table 5.1: Default values for TE, TR, TI for the different synthetic contrast weighted images.

Image Type	TE	TR	TI
T1W	10 ms	650 ms	
T2W	100 ms	4 500 ms	
PDW	10 ms	8 000 ms	
T1W FLAIR	10 ms	2 500 ms	1 050 ms
T2W FLAIR	90 ms	15 000 ms	3 100 ms
T2W STIR	100 ms	15 000 ms	300 ms
PDW STIR	5 ms	15 000 ms	300 ms
PSIR	10 ms	6 000 ms	500 ms
PSIR (vessel)	10 ms	8 000 ms	10 ms

Table 5.2: Default TE, TR, TI and TI2 values for DIR for the supported field strengths.

Image Type	Field Strength	TE	TR	TI	TI2
Double IR (WM supp)	1.5 T	100 ms	15 000 ms	420 ms	3 600 ms
Double IR (WM supp)	3.0 T	100 ms	15 000 ms	470 ms	3 750 ms
Double IR (GM supp)	1.5 T	10 ms	6 000 ms	700 ms	3 200 ms
Double IR (GM supp)	3.0 T	10 ms	6 000 ms	900 ms	3 600 ms
Double IR (fat supp)	1.5 T	5 ms	15 000 ms	280 ms	3 330 ms
Double IR (fat supp)	3.0 T	5 ms	15 000 ms	280 ms	3 330 ms

### 5.1.2 Adjusting the Image Contrast

SyMRI allows post-examination control of the contrast weighting in the synthetic contrast weighted images. This is done by adjusting TE, TR, TI and TI2. These options can be selected in the right-click menu, their toolbar icon or keyboard shortcuts as indicated below. When adjusting the contrast weighting by holding down the left mouse button, you can revert the changes by pressing «Esc» before releasing the left mouse button.

**Adjust TE** «Q» TE is adjusted by moving the mouse while holding down the left mouse button.

**Adjust TR** «W» TR is adjusted by moving the mouse while holding down the left mouse button.

**Adjust TR/TE** «E» TR and TE are adjusted simultaneously by moving the mouse while holding down the left mouse button.

**Adjust TI** «R» TI is adjusted by moving the mouse while holding down the left mouse button. The option is available when the inversion pre-pulse is activated.

**Adjust TI-TI** «T» TI and TI2 are adjusted by moving the mouse while holding down the left mouse button. The option is available when the double inversion pre-pulse is activated.

TE

TR

TE  
TR

TI

TI  
TI2

Note: When the Adjust T1 or Adjust T1-TI tool is selected, the mouse pointer will become a red circle when passing over contrast weighted images where the corresponding inversion pulse(s) are not active. Attempting to adjust T1 or T12 for such an image will trigger a dialog box, asking you if you want to activate the corresponding inversion pulse(s).

**ROI Right-Click>Use ROI to Null Tissue** The contrast weight of the image will be adjusted, making tissue within the ROI appear black in the image. This is done through an inversion pre-pulse with an optimized TI. To null tissue of a certain type: place and adjust the size of a ROI so that the ROI only includes tissue of the type you want to cancel out the signal from. This will adjust TI in the shown image. In a DIR image, the TI used to null CSF will be adjusted to keep nulling CSF, while the other TI will be optimized to null the tissue in the ROI. Functions to draw a ROI are described in Section 5.3.

Also, see Section 7.3.1 (Use ROI to Null Tissue).

**Right-Click>Image Settings>Inversion Prepulse** «Ctrl+I» An inversion pre-pulse will be activated in the image. The TI of the inversion pre-pulse is displayed to the left in the view and can be adjusted using *Adjust TI*.

**Right-Click>Image Settings>Double Inversion Prepulse** Activate double inversion pre-pulse in the image. The inversion times are displayed to the left in the view and can be adjusted using *Adjust TI-TI*.

**Right-Click>Image Settings>Enable PSIR** When PSIR is inactivated, negative values are visualized as the corresponding absolute value. Negative values occur in images with an inversion pulse, for example PSIR. In conventional MRI, only the absolute value can be displayed. Using SyMRI the negative values can be shown by activating PSIR. When PSIR is activated the string "PSIR (synthetic)" will be shown in the viewport HUD.

**Right-Click>Image Settings>Enter Time Parameters** Shows a dialog box with input fields for each time parameter in the current viewport. Clicking OK sets the values and closes the dialog box. Clicking Apply sets the values but does not close the dialog box. Clicking Cancel reverts back to the initial values before opening the dialog box and closes the dialog box.

### 5.1.3 Adjusting FLAIR Images

In a FLAIR image, fluid (CSF) is suppressed using a combination of an inversion pulse and an inversion delay such that the fluid signal is nulled. The combination of TR and TI used may vary across hospitals and it is important to be aware of the differences.

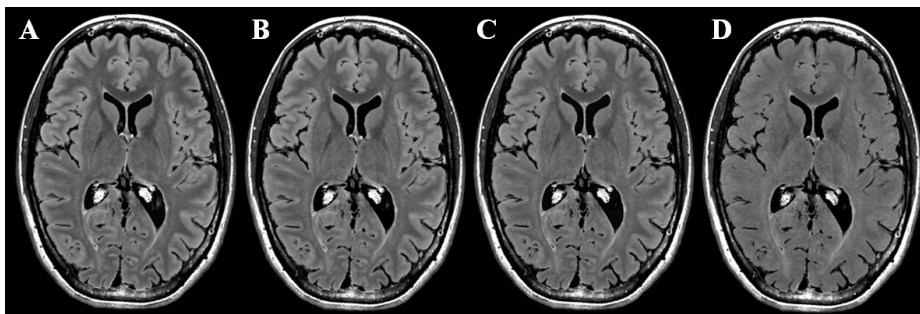


Figure 5.2: Various setting of T2W FLAIR on an MS patient (post-Gd).

A: TR/TI = 15 000/3 100 ms, B: TR/TI = 15 000/2 850 ms,  
C: TR/TI = 15 000/2 600 ms, D: TR/TI = 6 000/2 000 ms. The TE = 90 ms in all cases.

In Figure 5.2, a number of examples are given for optimization of FLAIR images: CSF is exactly nulled using the combination  $TR/TI = 15\ 000/2\ 850$  ms, resulting in very black CSF (Figure 5.2 B). However, a longer TI at the same TR still has a reasonably suppressed CSF but higher GM/WM contrast (Figure 5.2 A). A typical clinical FLAIR setting, using a short TR for scan time reduction, leads to degradation of GM/WM contrast but no change in scan time for synthetic MRI (Figure 5.2 D). SyMRI annotates an image as FLAIR when the TI value is near the value that suppresses CSF completely. For other TI values the image is annotated as IR. Also, see Section 7.3.1 (Use ROI to Null Tissue).

#### 5.1.4 Adjusting Double IR Images

In a Double IR image two tissues of choice are suppressed using a combination of two inversion pulses and inversion delay times. For example, WM, GM or fat can be suppressed simultaneously as CSF. In Figure 5.3a number of examples are provided together with the suggested settings.

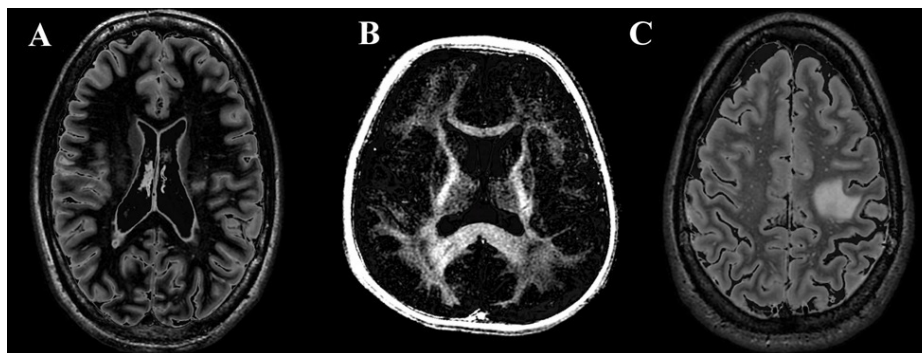


Figure 5.3: Example settings for Double IR at 3 T.

- A: A DIR where white matter and CSF are suppressed using  $TE/TR = 90/15\ 000$  ms,  $T11/T12 = 470/3\ 750$  ms (Female, 35 years),
- B: A DIR for pediatric applications, using  $TE/TR = 10/8\ 000$  ms,  $T11/T12 = 900/3\ 950$  ms (Male, 8 months) and
- C: A fat-suppression DIR similar to a T2W FLAIR image with fat-suppression, using  $TE/TR = 90/15\ 000$  ms,  $T11/T12 = 250/3\ 200$  ms (Female, 35 years).

To set two TI correctly is challenging and the “Use ROI to Null Tissue” (Section 5.1.2) function is recommended.

#### 5.1.5 Navigation Window for Contrast Images

The current contrast weighting of the image is displayed in the viewport in two ways; as a navigation window at the left side of the viewport and as text at the top of the viewport.

The navigation window (Figure 5.4) is displayed together with contrast weighted images. In the navigation window, the contrast weighting of the image is indicated by a yellow dot. When TR or TE is adjusted, use the navigation window to see how the changes affect the contrast weighting. It is also possible to adjust TE and TR by clicking on and dragging the yellow dot.

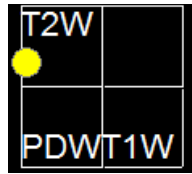


Figure 5.4: The circle in the navigation window indicates the contrast weighting of the image displayed in the viewport.

The text displayed at the top of the viewport changes depending on the current contrast weighting. For spin-echo images, it can be T1W, T2W, PDW or T1+T2 weight. If an inversion pre-pulse is active, the text will be IR, PSIR, STIR, T1W FLAIR, T2W FLAIR or FLAIR. If two inversion pre-pulses are activated, the contrast weighting will be DIR.

## 5.2 SyMaps

The quantification maps in SyMRI display the measured values for tissue properties per voxel in a dataset. The quantification maps for R1-relaxation rate, R2-relaxation rate and proton density (PD) are the foundation of synthetic MRI. In addition to these quantification maps the inverted relaxation times T1 and T2 are available.

### 5.2.1 Display of SyMaps

The quantification maps described below can be selected with the right click menu or the indicated keyboard shortcut.

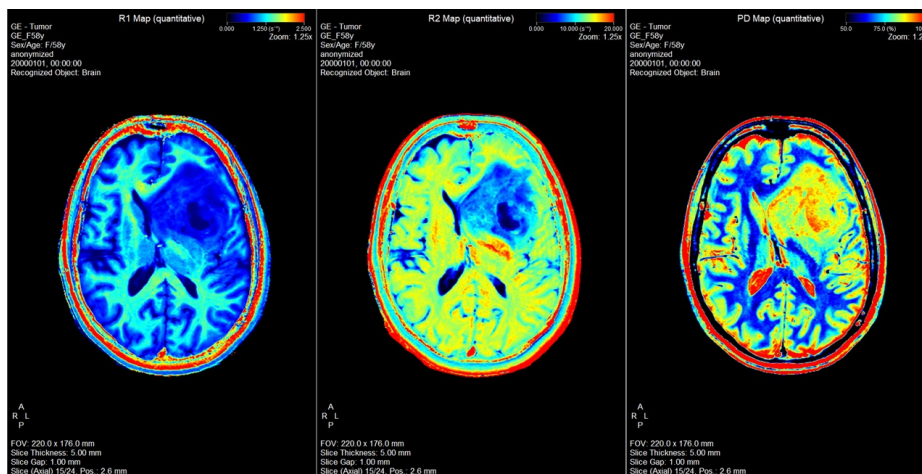


Figure 5.5: The SyMaps layout contains the SyMaps for R1, R2 and PD.

Each parameter represented as a SyMaps has a certain range within which it can be measured. This is known as the *dynamic range* of SyMRI. The dynamic range of each parameter is listed in Table 7.2.

**SyMaps Layout** «Ctrl+F» Display the default SyMaps layout which consists of 3 viewports with: R1-map, R2-map and PD-map.

**Right-Click>SyMaps>T1 Map** The T1 relaxation times are displayed in ms.

**Right-Click>SyMaps>T2 Map** The T2 relaxation times are displayed in ms.

**Right-Click>SyMaps>R1 Map** The R1 relaxation rate is displayed in  $s^{-1}$ ,

$$R1 = \frac{1}{T1}$$

**Right-Click>SyMaps>R2 Map** The R2 relaxation rate is displayed in  $s^{-1}$ ,

$$R2 = \frac{1}{T2}$$

**Right-Click>SyMaps>PD Map** The proton density is displayed in percentage units (pu) relative to a reference value of 100 pu, which is the proton density of water. Fat has higher proton density than water and is therefore displayed with values above 100 pu.



## 5.3 Region of Interest (ROI)

A ROI can be created to select an area in the image. Through the ROI, R1, R2, PD and the signal intensity within the marked region can be analyzed.

### 5.3.1 Create and Adjust ROI

Create a ROI by placing the mouse pointer on a selected spot, push down the left mouse button and draw desired shape and size, then release the mouse button. To activate the tool for drawing ROI, use any of the following alternatives:

**Freehand ROI** «F6» A ROI of free shape can be drawn. The ROI will be closed by a linear segment when the mouse button is released.

**Rectangle ROI** «F7» A rectangular ROI can be drawn. The rectangular ROI can be adjusted by marking one of the corners, holding the left mouse-button and adjusting the ROI to the desired size.



### 5.3.2 Information in ROI

Information related to the ROI is displayed in the upper right corner of the viewport (Figure 5.6), below the zoom and windowing. This block of text is referred to as *ROI Info*.

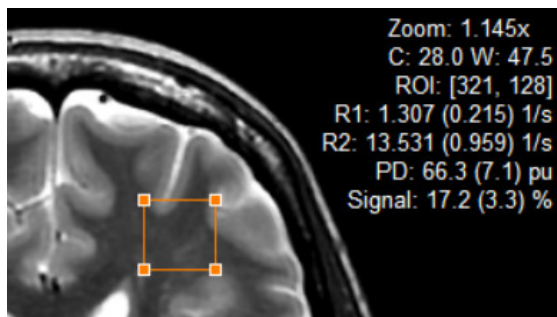


Figure 5.6: ROI and related information.

Two of the lines display the mean value and standard deviation for R1 and R2 by default. These two lines can be changed to display T1 and T2 instead.

The following information is displayed:

**Position:** The center coordinates of the ROI relative to the image are specified as "ROI [#,#]". For the free-hand ROI, the center is that of its bounding rectangle.

**R1 or T1:** Mean value and standard deviation of R1 or T1 within the ROI. For T1-maps, T1 is displayed and for R1-maps, R1 is displayed.

**R2 or T2:** Mean value and standard deviation of R2 or T2 within the ROI. For T2-maps, T2 is displayed and for R2-maps, R2 is displayed.

**PD:** Mean value and standard deviation of PD within the ROI.

**Signal:** Mean value and standard deviation of the voxel intensity. In synthetic images this will correspond to the signal strength of the simulated MRI signal.

Also, see sections 6.1, 6.4 and 6.7 for additional ROI operations available in SyMRI *NEURO*.

### 5.3.3 ROI Preferences

SyMRI supports several ROI per viewport and each ROI can be specific to a certain slice.

**Preferences > Multiple ROIs** Enables the user to have more than one ROI per viewport. Default is *off*. The active ROI has a solid outline while the other ROI have dotted outlines. The appearance of an active and inactive ROI can be seen in Figure 5.7.

When Multiple ROIs is *off*, a ROI is automatically deleted when a new ROI is drawn in the same viewport. When Multiple ROIs is *on*, existing ROI remain with a dotted outline.

**Preferences > ROI per Slice** Locks ROI to the specific slice where the ROI was drawn. Default is *off*.

When ROI per Slice is *off*, a ROI remains visible when scrolling through the slices. *ROI Info* is updated to display the statistics for the ROI in the currently visible slice. When ROI per Slice is *on*, a ROI is only visible in the slice where it was drawn. When scrolling to a different slice, the ROI will not be visible. *ROI Info* will only be visible in slices with a ROI.

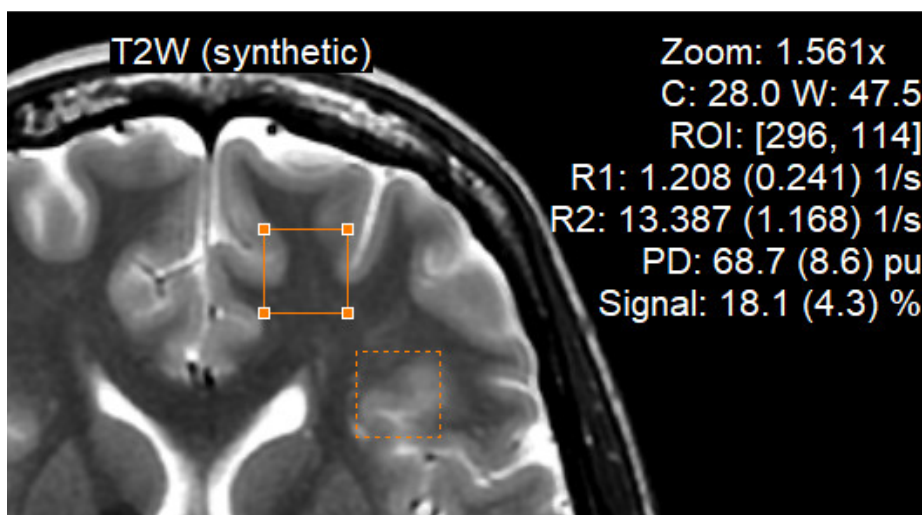


Figure 5.7: Two ROI in the same viewport. The upper ROI is active and therefore displayed with a solid line outline. The lower ROI is not active and therefore displayed with a dotted outline

Note: The information shown in the viewport is related to the ROI marked as active. A ROI is activated by clicking on it with the left mouse button. When a new ROI is drawn it is automatically marked as active. Which ROI is active may also be chosen by pressing «Tab» or «Shift+Tab».

### 5.3.4 Functions on a ROI

Right-clicking on a ROI will display the ROI right-click menu. This menu contains several operations for the ROI.

**ROI Right-Click > Cut ROI** «Ctrl+X» The ROI is copied and removed from the viewport.

**ROI Right-Click > Copy ROI** «Ctrl+C» The ROI is copied.



**ROI Right-Click>Paste ROI** «Ctrl+V» A previously copied ROI is pasted to the active viewport.

The key command «Ctrl+V» must be used to paste the ROI into a viewport without a ROI present. The key command «Ctrl+V» can also be used to paste the ROI as text. The text can then be copied and pasted into SyMRI to generate an identical ROI.

**ROI Right-Click>Disable ROI** «Backspace» The ROI is deleted from the viewport.

**ROI Right-Click>Show Plot** The plot will be hidden or shown.

**ROI Right-Click>Use ROI to Null Tissue** See Section 5.1.2 (Adjusting the Image Contrast). This function is automatically deactivated when changing slice.

## 5.4 Plots

With a ROI displayed in the image, a plot will automatically be displayed in the lower right part of the viewport. The values of the voxels in the ROI are plotted in the graph (Figure 5.8). The color indicates the amount of voxels that have a certain combination of values. Red indicates a relatively high number of voxels and blue a low number (Figure 5.9).

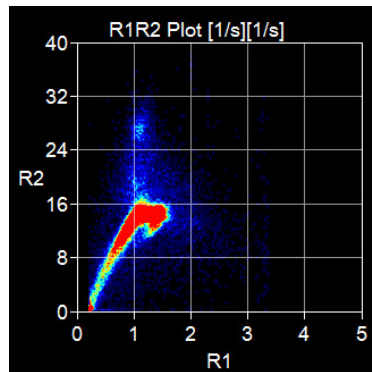


Figure 5.8: In the plot all voxels within the ROI are plotted.

The plot can be customized through the plot right-click menu. There are three types of graphs available in the plot right-click menu: R1-R2, R1-PD and R2-PD plots. The plot's axis and description change to match the selected plot type. There are two options available for the plot scale. One small plot scale which is adapted for values typically encountered in the brain. The other scale is a larger scale that provides a larger interval for the R1-R2, R1-PD and R2-PD axes. The small scale is the default and recommended scale for the brain anatomy, other anatomies will default to the large scale.

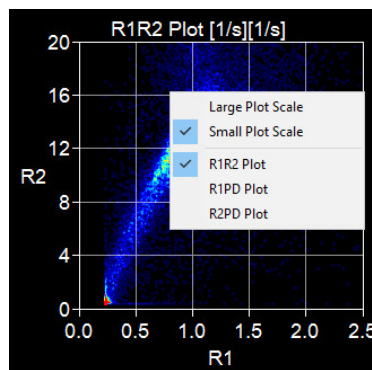


Figure 5.9: The plot right-click menu

## 5.5 Save and Load

Various types of images and image stacks can be saved from SyMRI to your workstation or PACS. Which save and load options are available depends on which deployment of SyMRI you are using.

**CAUTION:** Only the SyMaps and MDME source data can be re-opened by SyMRI. The saved contrast weighted image stacks cannot be opened by SyMRI, but can be opened by other DICOM image viewers.

---

**WARNING:** Observe that if the workstation is shut down or an error occurs while data is saved, the resulting data is not reliable. Remove any saved data and save again.



Note that for some deployments the save operation might only be partially completed when it has been completed in SyMRI. For instance, the Sectra PACS IDS7 deployment imports the saved data as a separate step.

### Plug-In

The Sectra PACS IDS7 starts the import when SyMRI initiates the save operation, however even after SyMRI exits the Sectra PACS IDS7 might still have images queued for import. If the Sectra PACS IDS7 crashes after SyMRI has exited you should inspect the saved data. Please refer to your Sectra PACS IDS7 documentation for details on how the Sectra PACS IDS7 import works.

---

### 5.5.1 Save SyMaps

The File-menu contains options for saving SyMaps.

Note: This function only available when using SyMRI with an MDME dataset.

Note: SyMaps are usually encrypted. If this is the case, there will be one file per slice. Encrypted SyMaps can be opened in a DICOM viewer, and will be displayed as a T2W series. The encrypted SyMaps can be opened correctly only in SyMRI.

MDME data for SyMRI typically contains over 300 images, while the quantified SyMaps contains one or three images per slice and approximately 20 to 300 images in total.

Never delete the original dataset! Observe that the original dataset might be needed to make reliable comparisons in future releases of the software.

### Plug-In

**File>Save SyMaps (T1T2PD) to PACS** A quantified image stack will be saved to PACS as a separate series in the current examination.

### StandAlone

**File>Save SyMaps (T1T2PD) to Local Folder...** A quantified image stack will be saved in DICOM-format to the selected folder on the local workstation. A **DICOMDIR** is created in the folder and the stack data is saved in a subfolder. If a **DICOMDIR** file already exists it is updated.

### Network

**File>Save SyMaps (T1T2PD) to PACS** A quantified image stack will be saved to the configured PACS via network as a separate series in the current examination.

### 5.5.1.1 Save Time

Before SyMRI is closed, saving the quantified SyMaps is recommended if you want to load the dataset at a later time. This is recommended because it is faster to load quantified SyMaps than starting from unquantified raw data (MDME).

### 5.5.2 Save a Screenshot of the Workspace

The File-menu contains options for saving a screenshot of the workspace.

Note: The screenshot image cannot be opened by SyMRI but can be viewed in other DICOM image viewers.

#### Plug-In

**File>Save Screenshot to PACS** Screenshot of the SyMRI workspace is saved to PACS as a new series in the current examination.

#### StandAlone

**File>Save Screenshot to File...** Screenshot of the SyMRI workspace is saved in DICOM format in the specified folder and given a specified file name on the workstation. A **DICOMDIR** is created in the folder and the file is saved in a subfolder. If a **DICOMDIR** file already exists, it is updated.

#### Network

**File>Save Screenshot to PACS** Screenshot of the SyMRI workspace is saved to the configured PACS via network, as a new series in the current examination.

### 5.5.3 Save Image Stacks

Note: The saved contrast weighted image stack cannot be opened by SyMRI, but can be viewed in other DICOM image viewers.

Note: The format of the saved images is different depending on the color setting. If color is *on*, the images will be saved as RGB-images. If color is *off* the images will be saved as grayscale. See Section 3.4.2.

#### 5.5.3.1 Save this Stack

The right-click menu for each viewport contains the option for saving the stack shown in the current viewport to a new series. The stack retains the settings from the viewport, e.g. windowing, contrast weighting and color settings.

#### Plug-In

**Right-Click>Save this Stack to PACS** Saves the stack to PACS as a new series in the current examination.

#### StandAlone

**Right-Click>Save this Stack to Local Folder...** Saves the stack in DICOM format to the specified directory on the workstation. A **DICOMDIR** is created in the folder and the stack data is saved in a subfolder. If a **DICOMDIR** file already exists it is updated.

Network

**Right-Click>Save this Stack to PACS** Saves the stack to the configured PACS via network, as a new series in the current examination.

### 5.5.3.2 Save all visible stacks

The toolbar and File-menu contains the option to save all currently visible viewports as new series. If raw data has been opened, SyMaps are also saved.



Plug-In

**File>Save all Visible Stacks to PACS** Saves all visible stacks (along with SyMaps if opened from raw data) to PACS as new series in the current examination.

StandAlone

**File>Save all Visible Stacks to Local Folder...** Saves all visible stacks (along with SyMaps if opened from raw data) in DICOM format to the specified directory on the workstation. A **DICOMDIR** is created in the folder and the stack data is saved in a subfolder. If a **DICOMDIR** file already exists it is updated.

Network

**File>Save all Visible Stacks to PACS** Saves all visible stacks (and SyMaps if opened from raw data) to the configured PACS via network, as new series in the current examination.

# Part III

SyMRI<sup>®</sup> *NEURO*

## 6 Functions in SyMRI *NEURO*

In SyMRI *NEURO*, the functionality to perform brain tissue characterization and brain volume measurement is available. Using the tissue properties R1, R2 and PD the tissue partial volumes are automatically measured. In this way, the intracranial volume (ICV) and the brain parenchymal volume (BPV), as well as white matter (WM), gray matter (GM), cerebrospinal fluid (CSF) and myelin correlated component (MyC) volumes are estimated. A user mask is available to measure volumes of user-defined selection of tissue.

**CAUTION:** The precision of these measurements depends on several parameters, including sequence settings, scanner manufacturer and field strength. Measuring the same patient again with the same sequence and the same MR scanner has a high repeatability. Any parameter change, however, such as changing the sequence or using another MR scanner, may decrease the precision. The brain volumes were validated at both 1.5 T and 3 T, but tolerance between scanners can always lead to differences. It is therefore recommended to use identical settings when doing longitudinal studies.

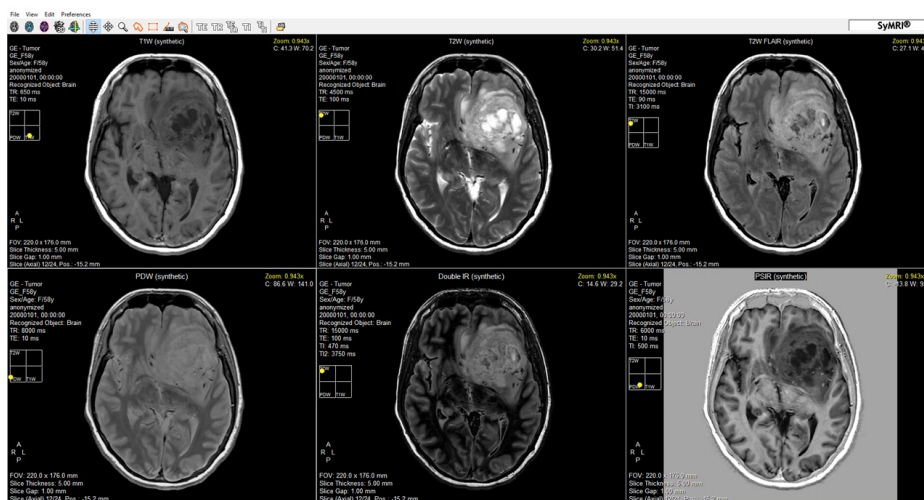


Figure 6.1: When SyMRI *NEURO* is started, six viewports are shown with synthetic contrast weighted images: T1W, T2W, T2W FLAIR, PDW, Double IR (WM supp), and PSIR. Use the keyboard accelerator «**Ctrl+R**» to return to this layout.

## 6.1 The Intracranial Mask

**WARNING:** Before analyzing tissue volumes, the user shall verify that the intracranial mask is anatomically correct. If the automatic generation is unsuccessful, the mask must be corrected by the user using the ROI functionality. Problems typically occur at the orbits and at the vertex and base of the skull.

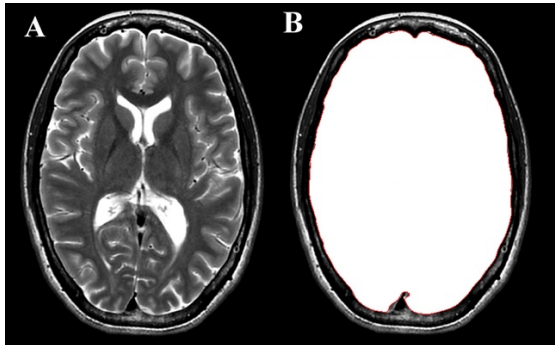


Figure 6.2: The intracranial mask defines the ICV in SyMRI *NEURO*. The red line indicates the outside of the ICV

The intracranial mask defines the intracranial volume in SyMRI *NEURO*. All segmentation and volume calculations are applied to the voxels included in the intracranial mask.

The intracranial mask is found by including WM, GM and CSF, followed by a region growing algorithm to ensure a contiguous volume. The edge of the mask is defined at PD = 50, based on the assumption that the ICV edge is halfway CSF (with PD = 100) and skull bone (PD = 0). Finally, the intracranial mask is made 0.5 mm smaller to remove the dura mater. Note that the red ICV edge line, as displayed in Figure 6.2, is plotted outside the ICV volume.

### 6.1.1 Adjusting the Intracranial Mask

The intracranial mask should be adjusted using the following steps:

1. **Show the intracranial mask** by selecting *Intracranial Mask* on the right-click menu, or by pressing «I».
2. **Scroll through all slices** and verify that the intracranial mask is correct.
3. **Add areas** to the mask by drawing a ROI around the area. Add the area to the mask by using *Add to Intracranial Mask* in the ROI right-click menu, or by pressing «Insert» on the keyboard. Parts of the ROI that extend outside the image are ignored when editing.
4. **Remove areas** from the mask by drawing a ROI around the area. Remove the area from the mask by using *Remove from Intracranial Mask* in the ROI right-click menu, or by pressing «Delete» on the keyboard. Parts of the ROI that extend outside the image matrix are ignored when editing.
5. **Save the intracranial mask** When the mask is correct, save changes using *Save Intracranial Mask* in the File-menu. When the mask is saved, the text below the segmentation table will be updated to: *"Intracranial mask last modified by 'name of user', 'date', 'time'"*.



#### Plug-In

The intracranial mask is not saved automatically when SyMaps are saved, but will be recalculated when the SyMaps are opened. It is possible to save the intracranial mask when SyMaps are opened. Once an intracranial mask has been saved with the SyMaps, this mask will be automatically loaded and displayed whenever the SyMaps are opened.

---

**WARNING:** The intracranial mask (Figure 6.2) corresponds to SyMRI *NEURO*'s definition of ICV. The tissue segmentation and volume estimation is performed for all voxels in the intracranial mask. Tissue outside the mask will not be included in the estimation of tissue volumes.

---



### 6.1.2 Additional Functions for the Intracranial Mask

Menu items for the intracranial mask are found in both the File- and Edit-menu:

**File>Save Intracranial Mask** The intracranial mask is saved and will be loaded automatically the next time the examination is used. This menu option is only available when a SyMaps dataset has been opened.

#### StandAlone

The information is saved to a presentation state DICOM-object file in a new series in the same folder the data-set was loaded from. When saving SyMaps the intracranial mask is saved to the saved SyMaps dataset.

#### Plug-In

The information is saved to a presentation state DICOM-object in PACS in a new series in the same examination as the loaded data-set.

#### Network

The information is saved to a presentation state DICOM-object in the configured PACS via the network in a new series in the same examination as the loaded data-set. When saving SyMaps the intracranial mask is saved to the saved SyMaps dataset.

**File>Load Intracranial Mask** The last saved intracranial mask is loaded.

**Edit>Recalculate Intracranial Mask** The intracranial mask is regenerated by the algorithm included in SyMRI *NEURO*.

**Edit>Clear Intracranial Mask** The intracranial mask is erased. All calculated volumes are set to zero. This function is useful when the user wants to define the brain from scratch.

**Right-Click>Show Intracranial Mask Edge** «H» Toggle show or hide the intracranial mask edge.

---

**WARNING:** To create accurate segmentation results the complete skull volume must be included in the intracranial mask. The intracranial mask shall not be used to estimate the volume of specific parts of the brain e.g. ventricles. To segment smaller regions of the brain the user mask (Chapter 6.4) can be used.

---



## 6.2 Tissue Segmentation

SyMRI automatically segments several tissue types. These include WM, GM, CSF and MyC. The automatic segmentation is done on all voxels contained within the ICV, defined by the intracranial mask.

The results of the tissue segmentations are presented in two ways, both as total volumes presented in the segmentation table (Section 6.5) and as segmentation maps (Section 6.3).

### 6.2.1 Brain Parenchyma and CSF Segmentation

Within the ICV, SyMRI segments tissue into four types, where each voxel within the ICV contains partial volumes in the interval 0 to 100 % of the voxel volume for each tissue type. The four tissue types in SyMRI's segmentation model are:

- **WM** , (Myelinated) White Matter
- **GM** , Gray Matter
- **CSF** , Cerebrospinal Fluid
- **Non-WM/GM/CSF**, tissue/fluid with characteristics which deviate from WM, GM and CSF

Note: Tissue segmented as Non-WM/GM/CSF should not be assumed to be pathological, it is merely the remainder of the voxel volume that was not segmented as WM, GM or CSF.

The tissue definitions of WM, GM and CSF were defined using analysis of R1, R2 and PD, as quantified by SyMRI. A neuroradiologist placed several ROI to mark areas of pure WM, GM and CSF in brain datasets of a number of adult healthy volunteers. Statistical analysis of the tissue properties within these ROI was used to define Gaussian distributions of R1-R2-PD for each tissue type, also referred to as tissue clusters.

---

**WARNING:** The tissue types are similar to but not identical to tissue according to histological classification. For example, in very young children, where white matter has not been fully myelinated, the tissue may be segmented as GM instead of WM. This happens since WM is modeled for *myelinated white matter* and the *unmyelinated* white matter more closely matches the tissue properties of GM.

---



---

**WARNING:** No comparison with biological white matter, gray matter or CSF has been made.

---



### 6.2.2 Myelin Correlated (MyC) Partial Volume Mapping

Separate from the segmentation maps for WM, GM, CSF and Non-WM/GM/CSF, SyMRI can also display a segmentation map for Myelin correlated component, MyC. This map has a high correlation to the partial volume of myelin in each acquisition voxel. Typical values for MyC are 0–8 % in cortical GM and 30–45 % in WM in adults. MyC is displayed in the image as a color scale from 0 to 40 %, but the values may be higher. Use the ROI function to verify MyC values in specific areas (see Section 6.7.1). In the human brain,

MyC typically does not exceed 50 %. MyC is not linearly proportional to the WM partial volume, e.g. MyC may vary in areas homogeneously segmented as 100 % WM.

MyC is based on another segmentation model than WM, GM, CSF and Non-WM/GM/CSF. The sum of the volumes of WM, GM, CSF and Non-WM/GM/CSF add up to the total ICV, without including MyC. Real myelin contains very thin layers of water, where relaxation occurs too fast to be directly observed by the SyMRI acquisition. Owing to magnetization exchange, however, the fast relaxing myelin water has an effect on the surrounding, slower relaxing water, comprising axonal, intracellular and extracellular water. Hence, the MyC values are based on the magnetic exchange between myelin and the surrounding water. When myelin is present, the observed relaxation rates of the surrounding water is higher, and the observed PD is lower, than a situation where no myelin is present. The MyC model contains four partial volumes, each with their own set of relaxation rates and proton density parameters. The four partial volumes are MyC, cellular partial volume, free water partial volume and excess parenchymal water partial volume. The sum of these four partial volumes is 100 % for each acquisition voxel. An algorithm is used to determine the optimal distribution of the four partial volumes to obtain the measured R1, R2 and PD values for each voxel in the parametric maps. The MyC partial volume has high R1 and R2 relaxation rates and low PD. The algorithm will find a higher contribution of MyC in WM than GM, since the relaxation rates in WM are higher than in GM, while the PD is lower in WM than in GM. Only the compartment containing MyC can be displayed in SyMRI. More details on the model for determining MyC are described in an open source journal which is freely downloadable via

<https://www.frontiersin.org/articles/10.3389/fneur.2016.00016/full>.

Left in figure 6.3 a synthetic T2-FLAIR of an MS patient is shown, where lesions are clearly visible. On the center image the MyC map is overlaid as color. A ROI function can show local MyC values, the title displays the volume per slice and in the segmentation table the entire volume is summed. On the right the MyC map is shown again, without background image and in B/W.

This image can be exported to PACS while keeping the numerical values.

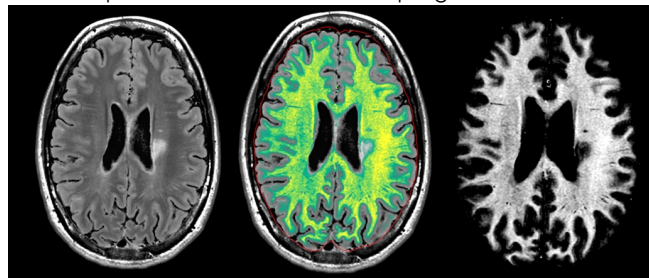


Figure 6.3: Example of MyC mapping.

Note: There is a high correlation between MyC values and optical density of histologically stained specimen with Luxol Fast Blue. Since the MyC represents the total magnetic effect of myelin on the R1, R2 and PD of surrounding tissue, however, it is not validated whether the partial volumes of MyC and real myelin are identical.

**CAUTION: Pathological changes in the brain tissue that causes R1, R2 and PD to deviate greatly from those found in non-pathological brain**

parenchyma may affect the MyC segmentation. One example of such a change is severe edema where the relaxation rates approach those of pure CSF.

## 6.3 Segmentation Maps

Segmentation maps are used to visualize the tissue content in each voxel within the ICV. As default, segmentation maps are displayed as overlays, colored images that are displayed over contrast weighted images. The contrast weighted image then acts as a frame of reference, while the colored overlay displays tissue content. The overlays have different color depending on which tissue they represent. Figure 6.4 displays five different segmentation maps, displayed as overlays over contrast weighted images. For WM/GM/CSF, the tissue maps model partial volume of healthy adult tissue, which means that the maps saturate at 100 % for tissue corresponding to pure WM/GM/CSF and (nonpathological) changes in e.g. myelination within pure white matter are not displayed since they are all mapped to 100 % in the WM map. The segmentation maps (Figure 6.4) are displayed on a background, which can be a contrast weighted image or solid black.

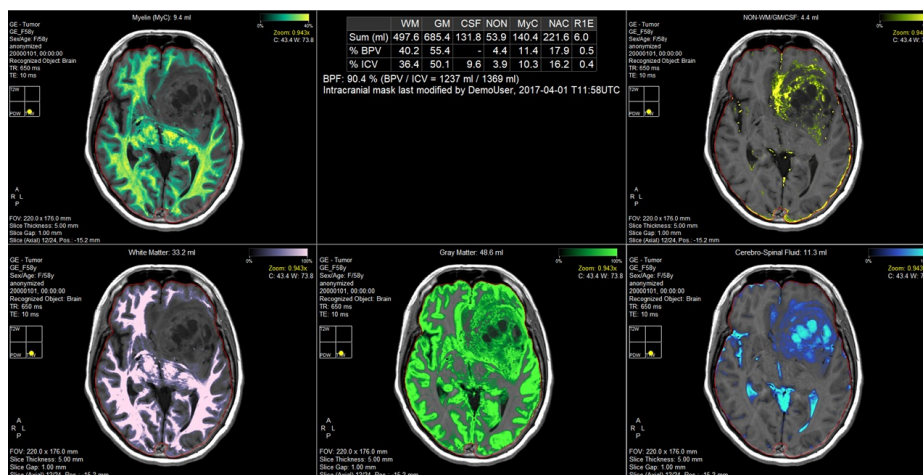


Figure 6.4: The segmentation layout consists of T1W images with MyC, WM, GM, CSF and Non-WM/GM/CSF as overlays together with the Segmentation Table

### 6.3.1 Display Segmentation Maps

The segmentation maps are activated in the Right-Click menu or through the associated keyboard accelerator.

**Segmentation Layout** «Ctrl+S» Display the default segmentation layout which consists of 6 viewports with: MyC, WM, GM, CSF and Non-WM/GM/CSF together with the Segmentation Table.

**Right-Click»White Matter** «7» The segmentation map for WM is displayed as an overlay.

**Right-Click»Gray Matter** «8» The segmentation map for GM is displayed as an overlay.

**Right-Click»Cerebro-Spinal Fluid** «9» The segmentation map for CSF is displayed as an overlay.

**Right-Click»Non-WM/GM/CSF** «0» The segmentation map for tissue that does not correspond to the definition of WM, GM or CSF is displayed as an overlay.



**Right-Click>Myelin (MyC) «Y»** The segmentation map for MyC is displayed as an overlay.

When a segmentation map is visible, the total volume of the segmented tissue within the current slice is shown in the top of viewport.

When a ROI is active and a segmentation overlay is visible, the volume within the ROI is shown in the upper right corner of the viewport. See section 6.7 for more information.

It is also possible to disable displaying the intracranial mask edge via menu option Show Intracranial Mask Edge or short key «H».

When a segmentation map is visible, a gradient colorbar appears in the upper right corner of the viewport. This bar shows how the color of the segmentation map relates to partial tissue volume. See figure 6.5 for an example of a colorbar.

### 6.3.2 Background Image for Segmentation

To view the segmentation maps on a linear scale from black to white, where white is 100 % and black is 0 % tissue content, the color can be turned off.

The color can only be turned off after the background is disabled.

When the background is disabled it is also possible to window the segmentation map in the same way as windowing a contrast image or a quantitative map (see section 3.5.3). The colorscale in the colorbar is updated when windowing.

**Right-Click>Disable Background «Ctrl+B»** Blank background is used. See figure 6.5.

**Right-Click>Image Settings>Color «C»** Disable color

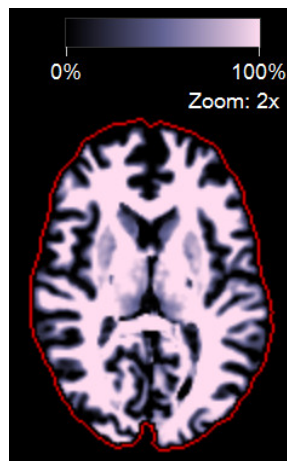


Figure 6.5: Example of segmentation overlay on blank background and a colorbar.

## 6.4 User Mask for User Defined Segmentation

The user mask gives the user the possibility to mark selected parts of the image stack and add it to a user mask. The user mask can include tissue both inside and outside the intracranial mask. The mask can be used to segment the volume of specific tissues, lesions and anatomies. The user mask is automatically displayed if SyMRI *NEURO* is started with a dataset that contains a previously saved user mask.

### 6.4.1 Display the User Mask

The user mask is activated in the right-click menu.

**Right-Click>Show User Mask** «U» The user mask is displayed in the viewport as an overlay.

### 6.4.2 Measuring Volume Using the User Mask

The main function of the user mask is to measure the volume of some user defined region(s). The volume is calculated using the voxel volume, including the underlying slice gap. The volume of the user mask can include all, or part of each voxel. The user mask volume in the current slice is displayed at the top of the view port(s) in which the user mask is displayed. Partial volume can be added to the user mask by copying partial tissue content from a segmentation map to the user mask (see 6.4.3). The partial volume of the segmented tissue is then added to the user mask volume. For example, if the partial volume of tissue in a voxel is 25 %, only 25 % of the voxel volume is added to the user mask volume.

User mask that only includes part of the voxel volume is partially transparent. More transparency represents lower partial volume. Options for color and background are identical to those for segmentation maps. See Section 6.3.2.

### 6.4.3 Create and Adjust User Mask

The user mask is created by adding or removing regions. This can be done by adding or removing regions marked by a ROI, by adding segmented tissue marked by a ROI, or by copying complete slices from a segmentation map.

#### 6.4.3.1 Edit User Mask Using a ROI

When using a ROI to edit the user mask, areas and tissues can be added or removed from the user mask via the ROI Right-Click menu. Start by placing a ROI around the area you wish to change.

**ROI Right-Click>Add to User Mask** All voxels within the ROI are added to the user mask. The total volume of each voxel is included in the user mask volume.

**ROI Right-Click>Remove from User Mask** Removes all voxels within the ROI from the user mask. The volume contributed by these voxels is removed from the user mask.

It is also possible to copy tissue from a segmentation map.

**ROI Right-Click>Copy Tissue in ROI to User Mask - Scale Up** All voxels containing more than 1 % of the tissue in the segmentation map is added to the user mask. The total volume of each voxel is included in the user mask volume.

**ROI Right-Click>Copy Tissue in ROI to User Mask - Partial Volume** Partial volumes are copied to the user mask.

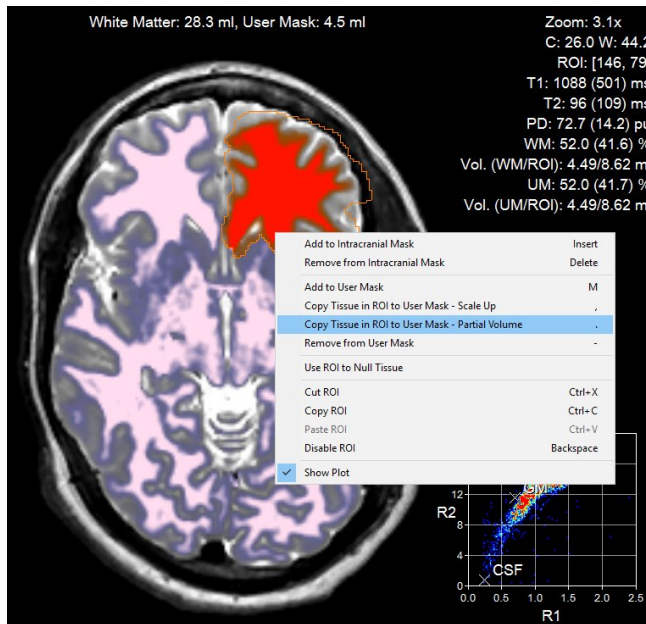


Figure 6.6:  
Copying with partial volume using a ROI with segmentation overlay

#### 6.4.3.2 Additional Ways to Adjust the User Mask

It is possible to add all segmented tissue, in one or more slices, to the user mask at once. This option appears on the right-click menu when a segmentation map is visible.

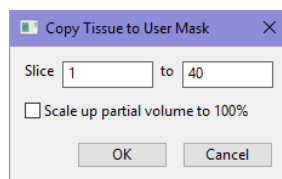


Figure 6.7: Complete slices of segmented tissue can be copied to the user mask.

**Right-Click>Copy Tissue to User Mask** Copy all tissue in the visible map, within a certain slice range, to the user mask. A dialog box appears, asking you to select the slice range (Figure 6.7). By default, partial tissue volume is added. Tick the box "Scale up partial volume to 100%" to include the total voxel volume of all voxels with more than 1 % partial tissue volume.

**Edit>Clear User Mask** All voxels are removed from the user mask and the user mask volume is set to zero.

#### 6.4.4 Save and Load the User Mask

It is possible to save the user mask, and to load a previously saved mask.



**File>Save User Mask** The user mask is saved. A dialog box appears asking you to give the user mask a name. It is possible to save several user masks with different names. The next time the data set is loaded in SyMRI *NEURO*, a list of all saved usermasks listed in the **DICOMDIR**-file will be displayed. If no other user mask is selected, or there is only one saved used mask, the last saved user mask will be displayed.

**StandAlone**

The information is saved to a presentation state DICOM-object file in a new series in the same folder the data-set was loaded from.

**Plug-In**

The information is saved to a presentation state DICOM-object in PACS in a new series in the same examination as the loaded data-set.

**Network**

The information is saved to a presentation state DICOM-object in the configured PACS via the network in a new series in the same examination as the loaded data-set.

**File>Load User Mask** The last saved user mask is loaded. Any current changes will be replaced by the last saved mask. The last saved user mask for a dataset will be loaded automatically when SyMRI *NEURO* is started.

## 6.5 Segmentation Table

The segmentation table provides volumetric information of all segmented tissue types in SyMRI *NEURO*, as well as the user mask. The total volumes for the whole brain are always shown, and the table can be expanded to show tissue volumes per slice.

	WM	GM	CSF	NON	MyC
Sum (ml)	491.9	868.7	106.4	28.7	158.7
% BPV	35.4	62.5	-	2.1	11.4
% ICV	32.9	58.1	7.1	1.9	10.6

BPF: 92.9 % (BPV / ICV = 1389 ml / 1496 ml)  
Intracranial mask generated by SyMRI NEURO

Figure 6.8: In the segmentation table, calculated volumes of WM, GM, CSF, Non-WM/GM/CSF and MyC are shown. The relative volumes, as a percentage of brain volume and ICV are also displayed.

**Right-Click>Segmentation Table** «6» Display the segmentation table in the viewport.

In the segmentation table the following volumes are displayed:

- **Sum (ml):** Total volume of WM, GM, CSF, Non-WM/GM/CSF and MyC in the intracranial mask.
- **% BPV:** Volumes of WM, GM, Non-WM/GM/CSF and MyC as percentages of the total brain parenchyma volume. Note that the percentages of WM, GM and NON sum up to 100 % while the percentage of MyC is a percentage of BPV independent of the other tissue types.
- **% ICV:** Volumes of WM, GM, CSF, Non-WM/GM/CSF and MyC as percentages of the total intracranial volume. Note that the percentages of WM, GM, CSF and NON sum up to 100 % while the percentage of MyC is a percentage of ICV independent of the other tissue types.

- **Brain parenchymal fraction (BPF).** The Parenchymal fraction is the sum of the volumes for WM + GM + Non-WM/GM/CSF divided by the total volume of the intracranial mask  

$$\left( BPF = \frac{BPV}{ICV} = \frac{ICV - CSF}{ICV} = \frac{WM + GM + NON}{ICV} \right).$$

Note: The percentage of MyC is not included in the calculation of BPF since it is independent of the other tissue types that completely partitions the content of ICV.

The name of the user that created or saved the intracranial mask is displayed below the table together with time and date.

### 6.5.1 Displaying Volumes per Slice

The segmentation table can be expanded to show the tissue volumes for each individual slice.

**Segmentation table right-click>Show Slices** expands the segmentation table to show the tissue volumes in each slice.

- **Slice** Volume of WM, GM, CSF, Non-WM/GM/CSF and MyC per slice.

### 6.5.2 User Mask in the Segmentation Table

If the user mask contains at least one voxel, a new column appears in the segmentation table named "UM". This column presents the user mask volume in the same way as WM, GM, CSF, etc. Both the volume per slice, total volume, and volume as percentage of ICV and BPV are displayed.

### 6.5.3 Exporting the Segmentation Table

**File>Save Segmentation Table as text...** The table with the segmentation results is saved locally on the workstation as a text file. In the file, patient name, examination ID and the versions of SyMRI used are saved.

**File>Save Segmentation Table as Structured Report...** The same command as above, but the table is saved as a DICOM object using the structured report format. Please refer to the DICOM Conformance Statement for SyMRI for further (technical) details on this format.

**Segmentation table right-click>Copy Segmentation Table** «Ctrl+C» Copy the segmentation table as text and lets you paste it into documents, emails etc.

## 6.6 Reference curves

The reference curves can be used to relate the brain segmentation volumes of the patient to a reference value, established from healthy subjects.

The reference curves can be shown via the right-click menu, or using the Reference Curves Layout icon. The layout shows reference curves for BPF, MyCPF, ICV, BPV and CSF.

The patient's value is indicated with a yellow dot. In the title of the reference curve plot the value is displayed, together with the reference value at that age and the standard deviations (sd) from the mean. See figure 6.9.



**Right-Click>Reference Curves** The reference curves sub-menu used to change visible reference curve.

**Right-Click on Axis** Change scale for axis.

The reference curves are based on data from:

1. McAllister A, Leach J, West H, Jones B, Zhang B, Serai S. Quantitative Synthetic MRI in Children: Normative Intracranial Tissue Segmentation Values during Development. *AJNR Am J Neuroradiol.* 2017;38:2364-2372,
2. Vågberg M, Ambarki K, Lindqvist T, Birgander R, Svenningsson A. Brain parenchymal fraction in an age-stratified healthy population - determined by MRI using manual segmentation and three automated segmentation methods. *J Neuroradiol.* 2016;43:384-391,

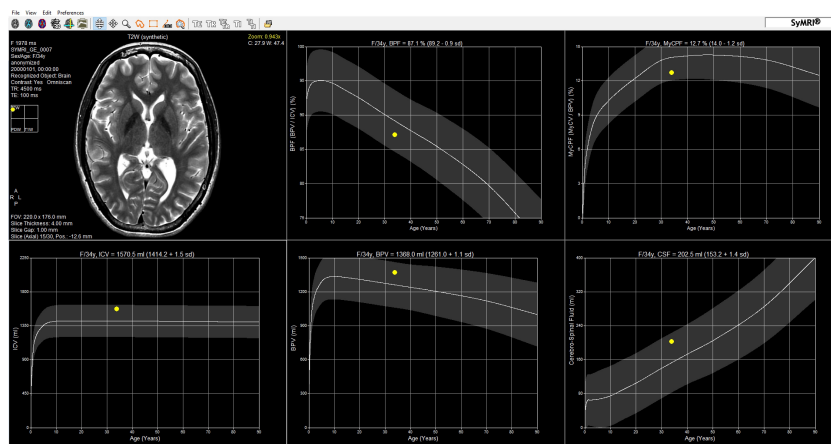


Figure 6.9: Screenshot of the reference curves layout.

## 6.7 ROI in SyMRI *NEURO*

In addition to the functionality that the ROI has in SyMRI *IMAGE* (Section 5.3) the ROI is used for three functions specific to SyMRI *NEURO*:

1. To adjust the intracranial mask (Section 6.1.1).
2. To adjust the user mask (Section 6.4.3).
3. Measuring the volume of segmented tissue.

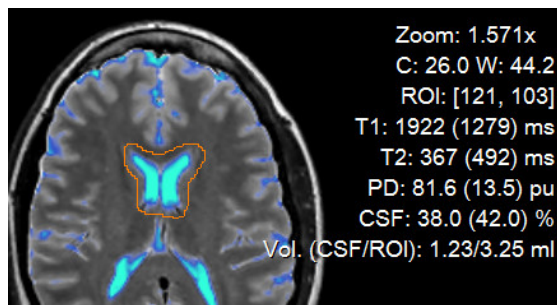


Figure 6.10: ROI and related information

### 6.7.1 Measuring the Volume of Segmented Tissue

When a segmentation map is visible, two new lines appear in the ROI Info text. The first line starts with the name of the visible segmentation map, e.g. WM or CSF. This line gives the mean and standard deviation of partial volume segmented tissue content within the ROI, as a percentage of the total ROI volume. The second line starts with Vol. and displays the volume of segmented tissue within the ROI, followed by the total volume of the ROI. The volumes are given in milliliters (ml).

If a user mask is visible, the volume information for the User Mask within the ROI will be displayed in the same manner in ROI Info. In case *both* a user mask and a segmentation overlay are visible, the volume information for the user mask is displayed below the volume information for the segmented tissue. See Figure 6.10.

Also note that the line displaying Signal is removed when the volume information for a user mask or segmented tissue is displayed.

## 6.8 Plots in SyMRI *NEURO*

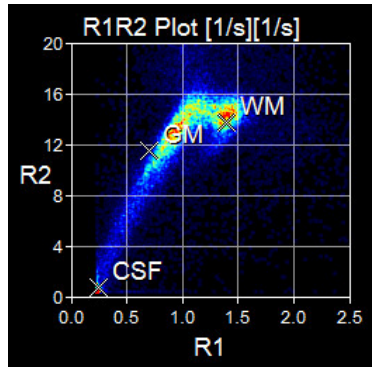


Figure 6.11: The ROI-plot in SyMRI *NEURO*.

The plots in SyMRI *NEURO* has the same functionality as in SyMRI *IMAGE* (Section 5.4) with the addition of displaying the mean value of the defined tissue types with a cross.

**Plot Right-Click>Brain Cluster** Reference clusters for CSF, WM and GM are displayed.

# Part IV

## Important Information

# 7 Warnings

Read all information in this chapter, it contains important information regarding safe usage of SyMRI.

## 7.1 Image Artifacts

The types of image artifacts appearing in conventional MR-images might also appear in synthetic MRI. Image artifacts such as motion artifacts, noise, folding etc. can appear for identical reasons. The effect on the images, however, might be different as the signal sampling is different.

### 7.1.1 Motion Artifacts

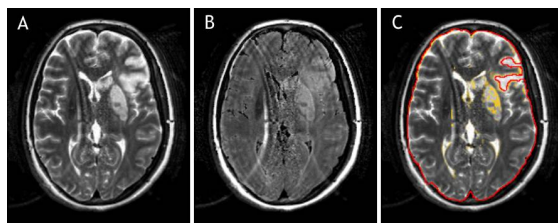


Figure 7.1: Patient motion might lead to folding and ghosting.  
A: T2W, B: FLAIR,  
C: T2W with Non-WM/GM/  
CSF-map

Just like conventional MRI, synthetic MRI is affected by movement. Movement results in artifacts that decrease the quality of the images. Synthetic images are made from several consecutive readings. Therefore, motion artifacts might appear different compared to conventional images. "Ghosting"-artifacts are similar to the ones seen in conventional images. In most cases the segmentation algorithm performs well, apart from the area affected by ghosting. Movement leads to fuzzy images with blurry edges.

---

**WARNING:** If the patient moves from her/his original position this will lead to errors in the transitions between tissues in the R2-map and severe measurement errors in the R1-map. This will result in unusable contrast images and poor segmentation. Figure 7.1 displays a stroke patient who has been moving during the examination.

---





## 7.2 Deviations - Synthetic MRI

The synthetic contrast weighted images are intended to look similar to conventional MRI images. There are however factors that make the synthetic images different from conventional ones.

### 7.2.1 Black Blood

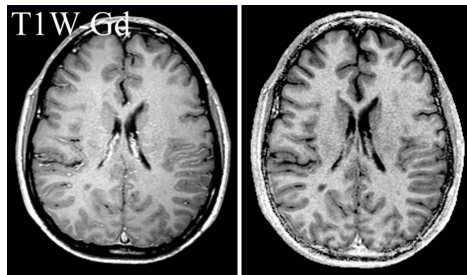


Figure 7.2: The signal from moving blood is suppressed in Synthetic images. For example arteries will not be enhanced in a synthetic T1W-image (to the right) after contrast agent has been administered, as they would in a conventional T1W-image (left).

The image acquisition for quantification suppresses the signal from moving blood. This results in synthesized images with "black blood"-features. Arteries will not be enhanced in the image, even if a contrast agent has been administered. Non-moving blood, caused by for example leakage, will enhance the T1-value in synthetic images. This also applies when contrast agent has accumulated within tissue.

---

**WARNING:** It is recommended to perform MDME acquisition prior to contrast agent administration.

---



### 7.2.2 Flow Sensitivity

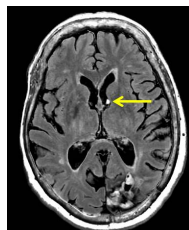


Figure 7.3: Hyper-intense signal due to a CSF flow jet.

Similar to conventional MRI, synthetic MRI may suffer from artifacts caused by flow. Examples are a band of ghosting artifacts due to a large artery or a flow jet in CSF (see Figure 7.3).

---

**WARNING:** Flow artifacts in particular are more pronounced for thinner slices. The thinner the slice, the more likely and pronounced the flow artifacts will be. Flow in arteries may also cause an apparent Hyper-intense Vessel Sign (also known as Arterial Hyperintensities) in FLAIR images.

---



### 7.2.3 Small Features Might Disappear

As in conventional MRI, small features will disappear if the resolution is too low for the particular feature of interest.

---

**WARNING:** In synthetic MRI, the calculations of R1, R2 and PD are based on a model assuming that the relaxation in each voxel is mono-exponential. Therefore, small structures not becoming the dominating component in any voxel might disappear and not be visible in the synthetic images.

---



## 7.3 Functions Specific for Synthetic MRI

### 7.3.1 Use ROI to Null Tissue

**WARNING:** Selecting large areas where there is more than one tissue type may cause the nulling algorithm to select a T1 which corresponds to neither tissue type. For instance, selecting a large ROI containing both WM and CSF in roughly equal proportions will cause incorrect nulling. See Figure 7.4 for an illustration of the effects. Best results are achieved by creating a ROI that contains only the tissue of interest.

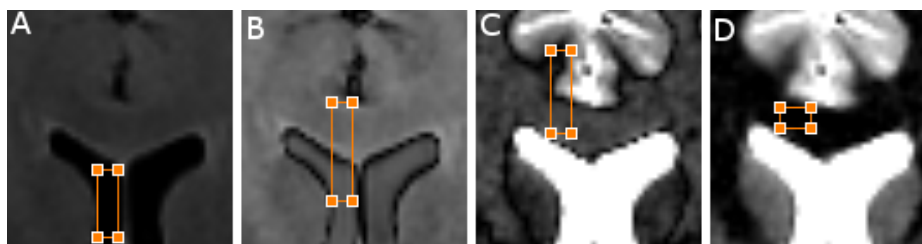


Figure 7.4: Effects of ROI placement for the null tissue function. Note that in 7.4d the caudate nucleus (*nucleus caudatus*) can be seen after the white matter has been cancelled.

### 7.3.2 Quantitative Measurements

The quantitative measurements of  $R_1$  (T1),  $R_2$  (T2) and PD have been verified in phantom studies. A set of phantoms with different combinations of  $R_1$  and  $R_2$  were scanned and quantified with SyMRI, and the results were compared to inversion recovery for  $R_1$  and multi-echo CPMG (Carr-Purcell-Meiboom-Gill) for  $R_2$ . For PD, a different phantom with various concentrations of heavy water for was used. A good correlation was shown for all parameters. The mean difference, standard deviation and verification range (interval) for the quantitative parameters, as well as different manufacturers and field strengths can be found in Table 7.1.

**Note:** Using SyMRI on MDME input data from Siemens 1.5 T scanners is not yet 510(k) cleared by the US FDA and thus currently not available for sale in the US unless limited to research or investigational use.

Table 7.1: Accuracy of SyMRI compared to gold standard on phantoms

Parameter	Interval	Philips		GE		Siemens	
		1.5 T	3.0 T	1.5 T	3.0 T	1.5 T	3.0 T
$R_1$	$0.5 - 3.33 \text{ s}^{-1}$	$3 \pm 6 \%$	$9 \pm 5 \%$	$5 \pm 6 \%$	$7 \pm 6 \%$	$-5 \pm 9 \%$	$-2 \pm 6 \%$
$R_2$	$2.5 - 25 \text{ s}^{-1}$	$11 \pm 2 \%$	$12 \pm 3 \%$	$11 \pm 8 \%$	$14 \pm 7 \%$	$13 \pm 4 \%$	$16 \pm 3 \%$
PD	20 - 100 pu	$8 \pm 9 \%$	$5 \pm 12 \%$	$3 \pm 6 \%$	$4 \pm 6 \%$	$6 \pm 11 \%$	$5 \pm 10 \%$

The results from the phantom studies are exemplified in Figure 7.5.

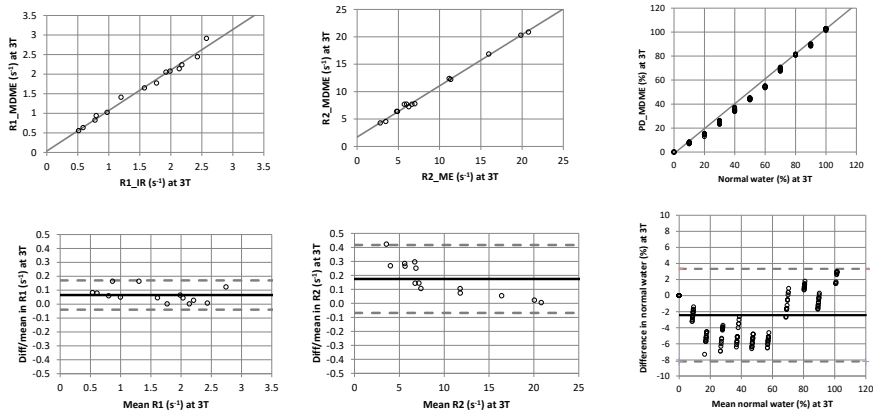


Figure 7.5: SyMRI quantitative measurements in comparison with gold standard quantitative measurements on phantoms for GE Discovery MR750 3 T.

### 7.3.2.1 Dynamic Range

The dynamic range of the measurements for the quantitative tissue parameters are listed in Table 7.2.

---

**WARNING:** Values outside the dynamic range will be truncated.

---



Table 7.2: Dynamic Range for SyMRI.

Value	Unit	Min	Max
T1	ms	250	4 300
R1	$s^{-1}$	0.23	4
T2	ms	10	2 000
R2	$s^{-1}$	0.5	100
PD	pu	0	160

### 7.3.3 Segmentation

The segmentation in SyMRI is based on pre-defined clusters (WM, GM, CSF) defined to represent white matter, gray matter and cerebro spinal fluid. Observe that pathological tissue can have values corresponding to any of the healthy definitions, in the same way healthy tissue might fall outside the defined tissue definitions. Partial volume is calculated through a simulation model (Bloch simulator) which aligns the values between two tissue clusters to estimated partial tissue type.

**CAUTION:** An absolute comparison between different brains, or a brain scanned at several points in time, must be made in comparable versions of SyMRI. Present Release Notes will state whether changes have been made compared to earlier versions of SyMRI.

---

**WARNING:** In voxels with partial volume of known tissue and partial volume of bone, air or other tissue with low PD, the segmentation may erroneously classify the voxel as Non-WM/GM/CSF instead of partial volume known tissue and partial volume Non-WM/GM/CSF.

---



---

**WARNING:** If there are mixtures of pathological tissues within a voxel the tissue within the voxel might be partially misclassified as a result of partial volume effects. For instance, a voxel containing two pathological tissues, where the mix intersects a known tissue, may cause misclassification as a partial volume of the known tissue.

---



---

**WARNING:** The administration of oxygen during MRI acquisition may increase R2 relaxation rate in brain tissue and hence affect the segmentation analysis.

---



### 7.3.4 Factors Affecting the Volume Estimations

Verification and validation have been performed to evaluate the repeatability of the volume measurements. The measurements are based on pre-defined tissue definitions.

#### 7.3.4.1 Accuracy

For volume estimations made in SyMRI, the volume for each separate voxel is calculated as the product of the voxel base area in the slice and the slice *spacing* (which includes the slice *gap*). The slice *spacing* is used instead of the slice *thickness* to include the complete brain volume. Figure 7.6 illustrates how the spacing, thickness and gap are connected.

**CAUTION:** SyMRI estimates the tissue content in the slice gap based on the two adjacent slices. The slice gap (Figure 7.6) shall be low, *but not less than 10 % of slice thickness*, to avoid large approximations and unacceptable cross-talk in the measurements. *The recommended slice gap is 25 % of the slice thickness*. See section 7.4.1.1 for information about the cross-talk and volume estimation trade-off.

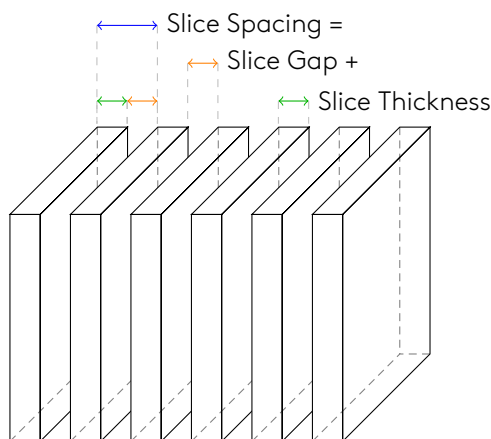


Figure 7.6: Illustration of the relationships between *Slice Spacing*, *Slice Gap*, and *Slice Thickness*.

### 7.3.4.2 Repeatability

To validate the repeatability of the segmentation method, in vivo measurements were performed on 7 healthy subjects (age 32–48 years) that were scanned using a GE 450W 1.5 T, a GE 750 3 T, a Philips Ingenia 1.5 T, a Philips Ingenia 3 T, a Siemens Aera 1.5 T, and a Siemens Prisma 3 T scanner. The subjects were scanned twice on each scanner, and were removed from the scanner between the acquisitions. Not all subjects were scanned using the Siemens scanners.

The repeatability of brain tissue volumes as a percentage of ICV, as well as the repeatability of BPF are reported in Table 7.3. The repeatability reported is the within-subject standard deviation, which means that the repeatability coefficient can be obtained by multiplying the values by 2.77.

Note: Using SyMRI on MDME input data from Siemens 1.5 T scanners is not yet 510(k) cleared by the US FDA and thus currently not available for sale in the US unless limited to research or investigational use.

Table 7.3: Repeatability for SyMRI.

	WM	GM	CSF	NON	MyC	BPF
GE 450W 1.5 T	0.46 %	0.58 %	0.26 %	0.13 %	0.14 %	0.26 %
GE 750 3 T	0.53 %	0.57 %	0.15 %	0.17 %	0.17 %	0.15 %
Philips Ingenia 1.5 T	0.51 %	0.73 %	0.12 %	0.19 %	0.20 %	0.12 %
Philips Ingenia 3 T	0.55 %	0.79 %	0.15 %	0.17 %	0.19 %	0.15 %
Siemens Aera 1.5 T	0.47 %	0.55 %	0.13 %	0.17 %	0.13 %	0.13 %
Siemens Prisma 3 T	0.38 %	0.35 %	0.15 %	0.12 %	0.13 %	0.15 %
Siemens Sola 1.5 T	0.36 %	0.41 %	0.14 %	0.11 %	0.15 %	0.14 %
Siemens Vida 3 T	0.44 %	0.46 %	0.10 %	0.19 %	0.11 %	0.10 %

**CAUTION:** To acquire comparable volume measurements, it is important that the same software version of SyMRI is used for all examinations. Volumetric measurements performed by different versions of SyMRI may be different due to changes in the software. Only compare datasets that have been quantified and segmented by the same version of SyMRI.

**CAUTION:** SyMaps from the same dataset but generated using different versions of SyMRI are not identical. Segmentation is performed base on the combination of R1, R2 and PD in each voxel. Therefore, any difference in the SyMaps may affect the segmentation and volume measurement. In short, using SyMaps generated by different versions of SyMRI may give different volumes even if they are both based on the same MDME acquisition data, and even if both are loaded into the same version of SyMRI.

## 7.4 Additional Warnings

### 7.4.1 System Requirements

---

**WARNING:** SyMRI is only to be used according to the system requirements. Any deviation from the system requirements may cause anything from very subtle changes, to rendering datasets unloadable by SyMRI. Attempting to load data that was not generated as specified by the system requirements violates the intended use of SyMRI.

---



#### 7.4.1.1 Slice Gap

The slice gap affects the quantification of R2 (T2) and volume estimations in segmentation.

**CAUTION:** The slice gap shall be *greater than 10 %* of the slice thickness to avoid unacceptable cross-talk between slices. Cross-talk affects the R2 (T2) quantification, and causes R2 to be overestimated. This will lead to an overestimation of MyC content.

*The recommended slice gap is 25 % of the slice thickness.* This gap is recommended since it provides a good trade-off between cross-talk and volume estimation error.

If the slice gap is wide, the estimation error for tissue volumes within the slice gap may increase. See Section 7.3.4.1.

#### 7.4.1.2 Slice Thickness

The slice thickness affects the quantification, resolution of small features and the risk of flow artifacts. The recommended slice thickness provides a good trade-off between resolution and flow artifacts.

If the slices are thin, flow artifacts become more likely. For more information, see Section 7.2.2.

If the slices are thick, the quantification will be affected by partial volume effects. For more information, see Section 7.2.3.

## 7.4.2 Follow-up Examinations Using Different Analysis Tools

---

**WARNING:** Different analysis tools use different models and methods of volume estimation, including different definitions of tissues. This results in systematic differences between different tools. Comparing the volumes given by SyMRI to other analysis tools in follow-ups of examinations can result in incorrect interpretation and/or conclusion. If comparing the segmentation volumes given by e.g. NeuroQuant, based on a T1W 3D acquisition, with the volumes given by SyMRI based on a 2D MDME acquisition, the volumes may differ for several reasons. Amongst others, the tissue definitions in NeuroQuant and SyMRI are different. E.g., SyMRI uses partial volumes, while NeuroQuant assigns one tissue type to each voxel.

In general, comparing data between different analysis tools should be avoided unless justification for the validity of comparing the different analysis tools can be provided.

SyntheticMR AB strongly recommend always using *the same analysis tool* to ensure that data from follow-up examinations can be compared, and to perform the comparisons using applicable clinical protocol and practice for follow-ups.

---

Note: SyMRI can only analyze tissue volumes on data from sequences that are approved for use with SyMRI. For a list of approved sequences with recommended sequence parameters please consult the System Requirement List.

## 7.4.3 Character Encoding

---

**WARNING:** If a dataset with non-west European characters for any of the textual information displayed in the HUD in SyMRI is loaded, this information will not be displayed correctly. For more information about what textual information is displayed in the HUD, see Section 3.2.1.1.

---

## 7.4.4 Anonymization Software

---

**WARNING:** SyMRI uses certain private DICOM tags to be able to load datasets (see DICOM conformance statement). If anonymizing software is used, check that these tags are preserved.

---





## 7.5 Scanner Specific Warnings

### 7.5.1 GE

The following warnings only applies to usage on GE MRI scanners.

#### 7.5.1.1 Anonymization

**CAUTION:** GE uses private DICOM tags to store a few acquisition parameters required by SyMRI. Use anonymization tools provided by GE to preserve these tags in anonymization. Datasets where the tags are missing can not be opened by SyMRI.

The tags required by SyMRI are listed in the DICOM conformance statement. If an anonymization tool not provided by GE is used, please confirm that the required tags are preserved.

### 7.5.2 Philips

The following warnings only applies to usage on Philips MRI scanners.

#### 7.5.2.1 Multiple Packages

Note that the System Requirement List requires that the scanner does not generate multiple packages.

On Philips MRI scanners with Philips software version 4.x.y (only available for 3T) or 5.1.x (available both for 1.5T and 3T) the scanner console will show that it is going to generate multiple packages but it does not store this information in the DICOM-objects; this means that it is not possible for post-processing software, such as SyMRI, to detect if the scanner has generated multiple packages after the scanning. Therefore, it is up to the scanner operator to ensure that the scanner reports that it will not generate multiple packages.

From Philips software version 5.2.x multiple packages will generate a conflict in the exam card and scanning can not be started.

---

**WARNING:** Loading multiple package data might cause unpredictable results which might not be apparently incorrect.

---



#### 7.5.2.2 Slice Order

**CAUTION:** On Philips MRI scanners with Philips software version 4.x.y (only available for 3T) or 5.1.x (available both for 1.5T and 3T) make sure that slice order in the exam card is set to Ascending or Descending and not Default or Interleaved.

The latter will cause quantification of T1 to oscillate between slices which can be seen by selecting a large ROI and scrolling through the slices. The curve for normal brain tissue in R1-R2-plot will then shift to the right/left almost every other slice. From Philips software version 5.2.x the slice order options Default and Interleaved have been removed.

### 7.5.3 Siemens

The following warnings only applies to usage on Siemens MRI scanners.

**Note:** Using SyMRI on MDME input data from Siemens 1.5 T scanners is not yet 510(k) cleared by the US FDA and thus currently not available for sale in the US unless limited to research or investigational use.

#### 7.5.3.1 Patient Position

**CAUTION:** Always position the patient so that the imaged anatomy is placed at the ISO-center of the scanner. If the Patient is incorrectly placed, the B1-field will be sub-optimal and the quantification may not be performed correctly.

MRI scanners have optimal imaging conditions at the center of the coil, the ISO-center of the magnetic field. Using the laser to position the patient correctly will ensure optimal quantification, and thereby optimal image quality in the synthetic images.

# Part V

## Appendix

## 8 Keyboard Shortcuts

A listing of all available keyboard shortcuts and their function available in SyMRI.

### macOS

On macOS the following shortcut keys are different from Windows:

- ⌘ instead of **Ctrl**.
- F13 instead of **Insert**.

Keyboard shortcut	Function
1	T1W
2	T2W
3	T2W FLAIR
4	STIR
5	PSIR
Page Up	Go to preceding slice
Up-arrow	Go to preceding slice <sup>1</sup>
Page Down	Go to succeeding slice
Down-arrow	Go to succeeding slice <sup>1</sup>
Home	Go to first slice
End	Go to last slice
Ctrl+R	Contrast Layout
Ctrl+F	SyMaps Layout
A	Auto Scale
Ctrl+Shift+I	Invert Grayscale
C	Color
Ctrl+Z	Zoom in
Ctrl+Shift+Z	Zoom out
Ctrl+P	Bicubic interpolation
Ctrl+X	Cut ROI
Ctrl+C	Copy ROI
Ctrl+V	Paste ROI

<sup>1</sup> This shortcut is not available in SyMRI Plug-In since the Sectra PACS IDS7 host window uses it.

Backspace	Disable ROI
Tab	Activate next ROI <sup>2</sup>
Shift+Tab	Activate previous ROI <sup>2</sup>
Ctrl+1	1 Viewport
Ctrl+2	2 Viewports
Ctrl+3	3 Viewports
Ctrl+4	4 Viewports
Ctrl+5	5 Viewports
Ctrl+6	6 Viewports
F4	Paging Tool
F5	Pan Tool <sup>1</sup>
P	Pan Tool
F6	Freehand ROI
F7	Rectangle ROI
Z	Zoom Tool
Q	Adjust TE
W	Adjust TR
R	Adjust TI
T	Adjust TI-TI <sup>3</sup>
X	AutoROI Tool <sup>4</sup>
Ctrl+I	Inversion Prepulse
Ctrl+0	Open...
Ctrl+Shift+S	Save all Visible Stacks to Local Folder...

**Pan Tool** is also available by holding the «Shift» key pressed.

**Zoom Tool** is also available by holding the «Ctrl» key pressed.

«S» and «F» can be used together with **Adjust TR**, **Adjust TI** or **Adjust TI-TI** to decrease or increase the rate by which the time parameters are changed when the cursor is moved. They can also be used to scroll faster between slices and when adjusting the windowing.

When the *User Mask* feature is available the following keyboard shortcuts can be used.

Keyboard shortcut	Function
U	Show User Mask

<sup>2</sup> For Multiple ROIs

<sup>3</sup> *DIR*-license is required.

<sup>4</sup> *AutoROI Tool*-license is required.

Ctrl+B                      Disable Background

With SyMRI *NEURO* the following keyboard shortcuts are also available.

<b>Keyboard shortcut</b>	<b>Function</b>
Ctrl+S	Segmentation Layout
Ctrl+C	Copy Segmentation Table
6	Segmentation Table
7	White Matter
8	Gray Matter
9	Cerebro-Spinal Fluid
0	NON-WM/GM/CSF
Y	Myelin (MyC)
I	Intracranial Mask
H	Show Intracranial Mask Edge
Ctrl+B	Disable Background
Insert	Add to Intracranial Mask
Delete	Remove from Intracranial Mask
M	Add to User Mask
Comma (,)	Copy Tissue in ROI to User Mask - Scale Up
Period (.)	Copy Tissue in ROI to User Mask - Partial Volume
Subtract (-)	Remove from User Mask

# SyMRI®

SyntheticMR AB

2019-12-16

## SyMRI® 11

Version 11.2.2

UDI: (01)07340024700030(8012)11.2.2

SyntheticMR AB

Storgatan 11

SE-582 23 Linköping, SWEDEN

[www.syntheticmr.com](http://www.syntheticmr.com)

[support@syntheticmr.com](mailto:support@syntheticmr.com)